

Rezūm therapy for ≥ 80 -mL benign prostatic enlargement: a large, multicentre cohort study

Dean Elterman, Naeem Bhojani, Christopher Vannabouathong, Bilal Chughtai, Kevin C. Zorn

Published April 24, 2022, in BJUI International

Limitations of this report include:

- Some patients had missing outcome data over time, requiring complete case analysis, which may limit statistical power or introduce bias.
- Treatment allocation was not blinded.
- Outcomes were measured over 12 months; longer-term durability beyond a year remains unknown.
- The study was conducted at two high-volume centers with experienced operators, which may not reflect outcomes in community or lower-volume settings.

To evaluate the efficacy and safety of Rezūm therapy in benign prostatic hyperplasia (BPH) patients with prostates ≥ 80 mL, a prospective registry was established at two high-volume Canadian centres.

A cohort of 83 patients with a large prostate (median prostate volume 100 mL) were included and followed for 12 months with assessments of urinary function, symptoms, and sexual health.

- IPSS score decreased by 24%, 57%, and 59% at 1 month, 3 months, 12 months, respectively ($p < 0.001$).
- IPSS-QoL decreased by 27% and 70% at 1 month and 12 months, respectively ($p < 0.001$).
- Qmax increased by 38%, 55%, and 59% at 1 month, 3 months, 12 months, respectively.
- Post-void residual decreased by 62% at 12 months.
- IIEF-15 increased by 15% at 6 months.
- Ejaculatory function improved initially, with 22.4% gain at 1 month, and only 3.6% of patients reported reduced or absent ejaculation.
- Re-catheterization was needed in 12 patients (14.5%) after an initial failed voiding trial, but all patients were voiding spontaneously within a month.



Rezūm™ System Brief Summary

Scan the QR code, visit
bostonscientific.com/rezum-risks

Disclaimers

Results from case studies are not necessarily predictive of results in other cases. Results in other cases may vary. The content of this article/publication is under the sole responsibility of its author/publisher and does not represent the opinion of Boston Scientific.

Please be aware that this summary has been created by Boston Scientific, and modifications were made from the original publication.


Rx Only.

CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labelling supplied with each device. Information for use only in countries with applicable health authority registrations. Material not intended for use in France.

All trademarks are the property of their respective owners.

© 2026 Boston Scientific Corporation or its affiliates. All rights reserved. URO-2215917-AB JUN 2025

CE 0344

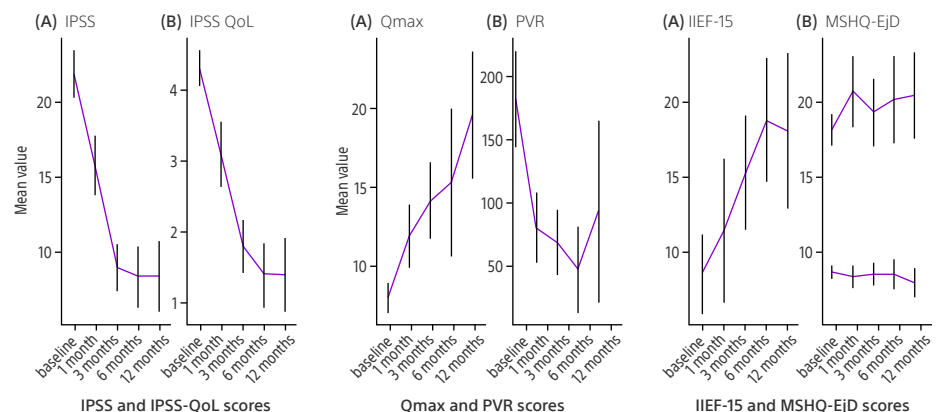

This is the first prospective, real-world study demonstrating that Rezūm therapy is safe and effective in large prostates (≥ 80 mL). Improvements in the above outcomes occurred as early as 1 month and up to 12 months after surgery. Overall, Rezūm therapy was shown to be a safe and effective procedure for patients with large prostate glands.

Full Article Reference: [Rezūm therapy for \$\geq 80\$ -mL benign prostatic enlargement: a large, multicentre cohort study](#)

[License Link](#)