

# Hurricane™ RX

Biliary Balloon Dilatation Catheter

**Boston  
Scientific**  
Advancing science for life™



*RX Biliary System™*

**We're in!™**

# Hurricane™ RX

## Biliary Balloon Dilatation Catheter

The Hurricane RX Biliary Balloon Catheter is a double-lumen catheter for endoscopic dilation of strictures of the biliary tree and the Sphincter of Oddi. The balloon is available in a range of sizes to accommodate patient anatomy and stricture size.

*RX Biliary System™*

# We're in!™


The RX Biliary System helps put the physician in control of the guidewire. Designed to facilitate access and control during ERCP, the system enables the physician to secure the guidewire to maintain access during exchanges — particularly helpful when working in tortuous anatomy — and provides a tactile feel while manipulating the guidewire.

Since its introduction, the RX Biliary System continues to provide utility with the development of new technologies and the continual improvement of devices for diagnostic and therapeutic ERCP.

### Excellent Stricture Passage

Low-profile, 5.8Fr balloon shaft is designed to reduce resistance and increase pushability when crossing strictures

- 4Fr tip is designed for less traumatic advancement across the stricture
- Stainless steel stiffening stylet inserted in the proximal portion of the catheter shaft to augment pushability
- RX catheter design allows guidewire locking to secure access

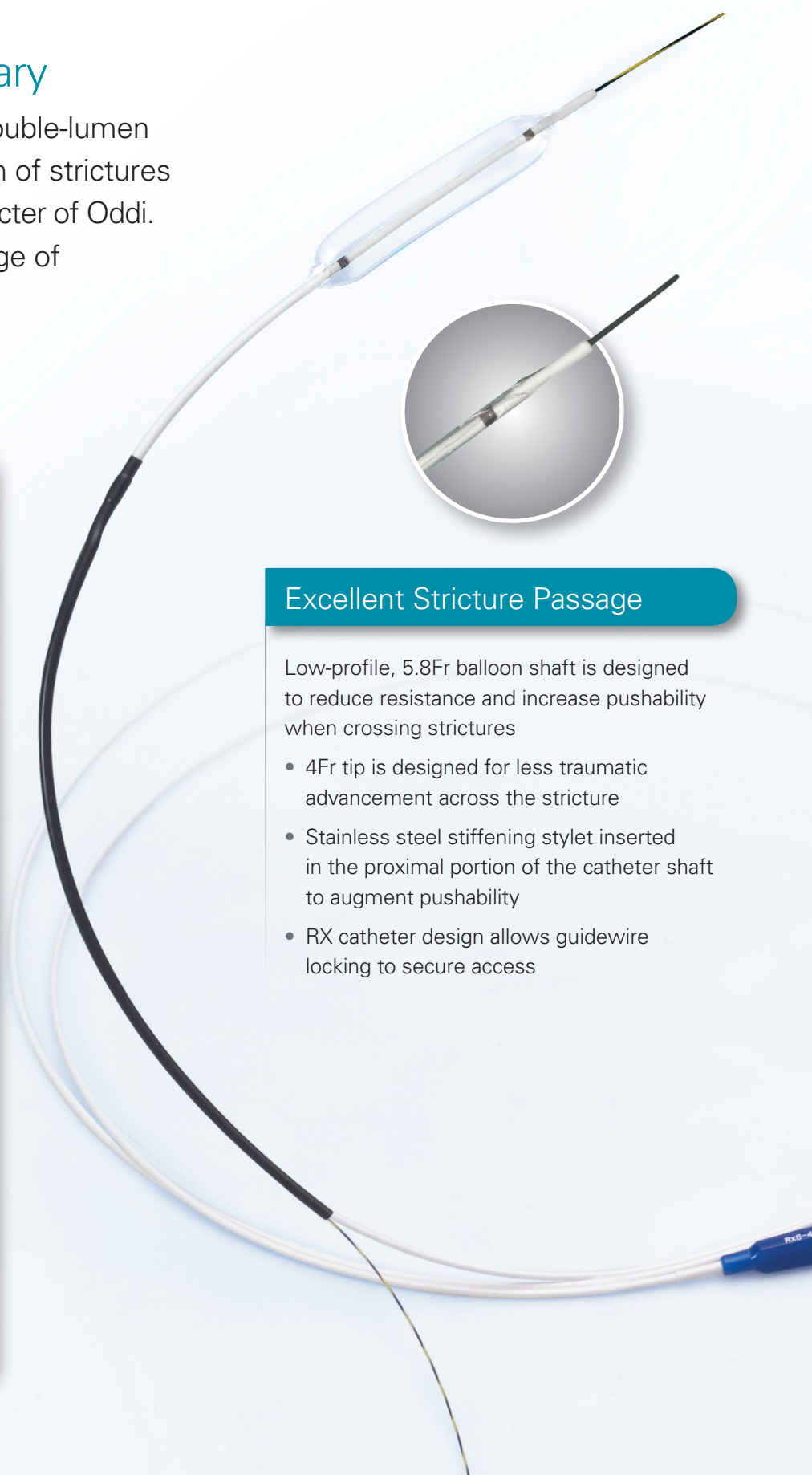

## Designed for Accurate Positioning


Hurricane™ RX Balloon offers multiple options for visualization

- Two radiopaque markers aid in balloon positioning
- Contrast may be injected with guidewire in place

## Range of Balloon Sizes

Options of balloon diameters and lengths facilitate selection of appropriately sized balloon for patient anatomy and stricture size

- 2cm and 4cm lengths available
- 4mm, 6mm, 8mm and 10mm balloon diameters available

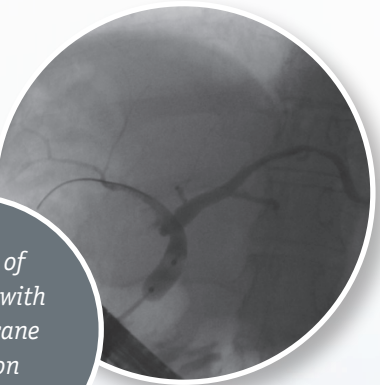

*Dilation of a stricture with the Hurricane Dilatation Balloon*

*Image courtesy of Mario Traina, MD, ISMETT Hospital, Palermo, Italy*

## Balloon Achieves Pressure Up to 12Atm

- Long 4cm balloon length is designed to offer stability during stricture dilation
- Balloon is recommended for use with the Encore™ 26 Inflation Device which has been developed for rapid deflation to potentially reduce patient discomfort



*Image courtesy of Sri Komanduri, MD, Feinberg School of Medicine, Northwestern University*

# Hurricane™ RX

## Biliary Balloon Dilatation Catheter



### Hurricane RX Biliary Balloon Dilatation Catheter

Order Number	GTIN	Balloon Inflated O.D. (mm)	Balloon Length (cm)	Catheter O.D. (Fr)	Catheter Length (cm)	Maximum Inflation Pressure (ATM)	Minimum Channel Size (mm)	Recommended Guidewire (in)	Packaging
M00545890	08714729283805	4	2	5.8	180	12	3.2	.035	Each
M00545900	08714729283812	4	4	5.8	180	12	3.2	.035	Each
M00545910	08714729283829	6	2	5.8	180	11	3.2	.035	Each
M00545920	08714729283836	6	4	5.8	180	11	3.2	.035	Each
M00545930	08714729283843	8	2	5.8	180	11	3.2	.035	Each
M00545940	08714729283850	8	4	5.8	180	11	3.2	.035	Each
M00545950	08714729283867	10	2	5.8	180	8	3.2	.035	Each
M00545960	08714729283874	10	4	5.8	180	8	3.2	.035	Each

\*Recommended Guidewire: .035" Hydra Jagwires Guidewire.

### Encore™ 26 Inflation Device

Order Number	GTIN	Description	Packaging
M00566670	08714729755241	Encore 26 Inflation Device	Each

**Boston Scientific**

Advancing science for life™

Boston Scientific Corporation  
300 Boston Scientific Way  
Marlborough, MA 01752-1234  
www.bostonscientific.com

**Ordering Information**  
**1.888.272.1001**

©2015 Boston Scientific Corporation  
or its affiliates. All rights reserved.

ENDO-334301-AB November 2015

All trademarks are the property of their respective owners.

Indications, Contraindications, Warnings and Instructions for Use can be found in the product labeling supplied with each device.

**Caution:** Federal (U.S.) law restricts this device to sale by or on the order of a physician.