


ATHLETIS™ Ultra-High-Pressure Balloon

ULTRA-HIGH PRESSURE. ULTRA-LOW PROFILE.


Optimized Deliverability to Reach and Treat Challenging Lesions

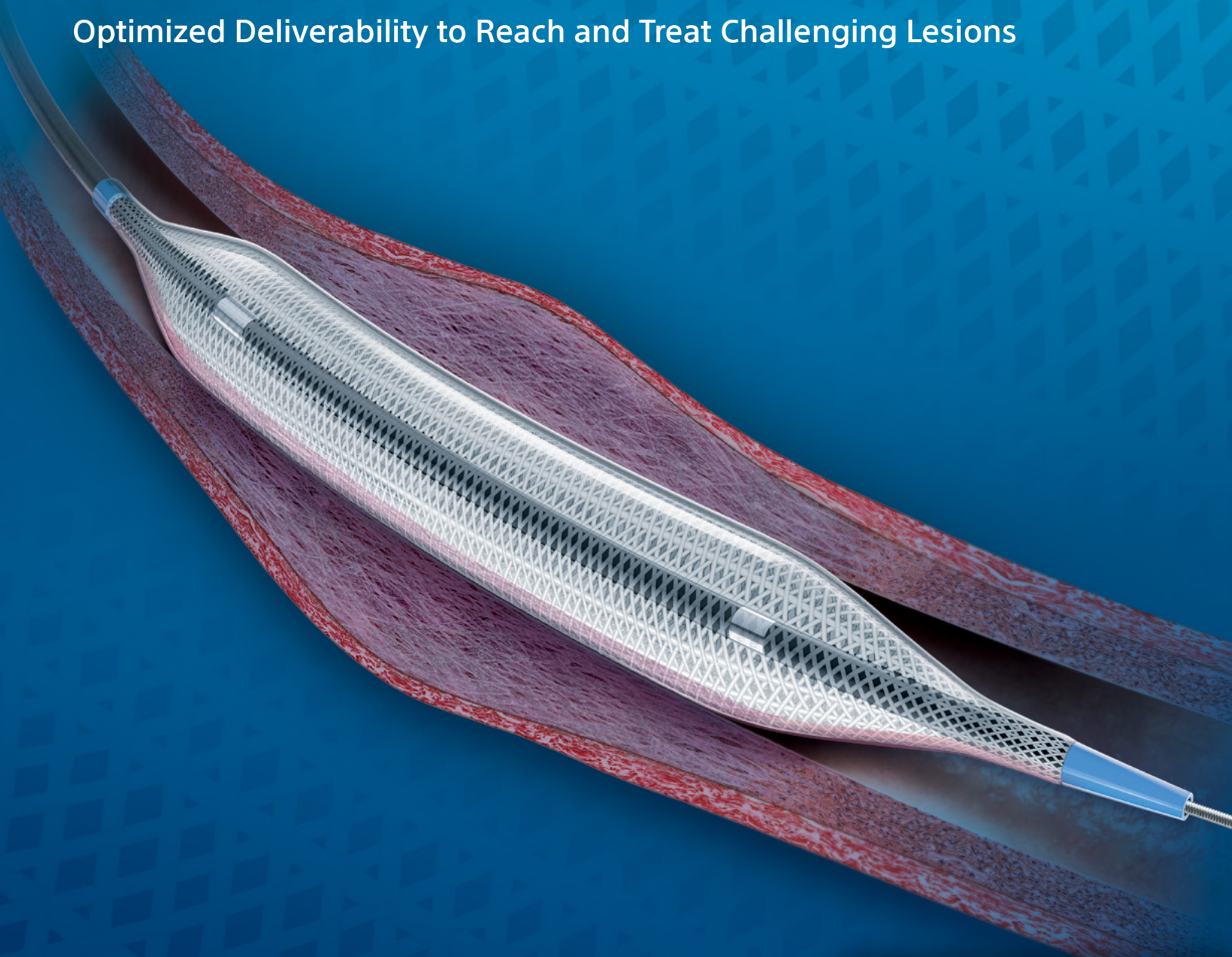

ULTRA-HIGH PRESSURE

Athletis Ultra-High-Pressure Balloon's proprietary braided balloon design enables optimized deliverability and ultra-high pressures up to 40 ATM. The unique, non-compliant balloon design enables Athletis Balloon Catheter to open to the stated diameter even in highly resistant lesions.

The braiding matrix was designed using a liquid-crystal polymer material that was used by NASA in the Mars Rover parachute.

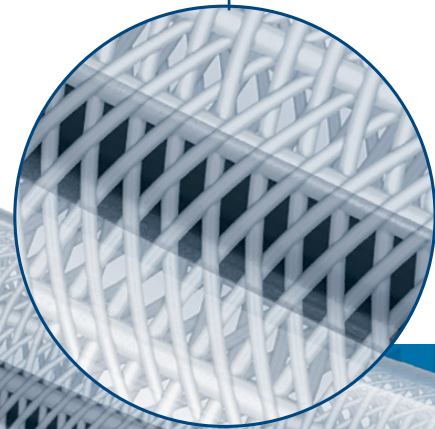


OPTIMIZED DELIVERABILITY

With an ideal blend of track and a low crossing profile, Athletis Ultra-High-Pressure Balloon is designed to reach highly calcified, long, stenotic, or resistant fibrotic lesions, including those commonly found in AV access procedures.

ULTRA-LOW PROFILE

Athletis Ultra-High-Pressure Balloon features a class-leading low lesion entry profile with an 0.042" flexible, tapered tip over an 0.035" guidewire.

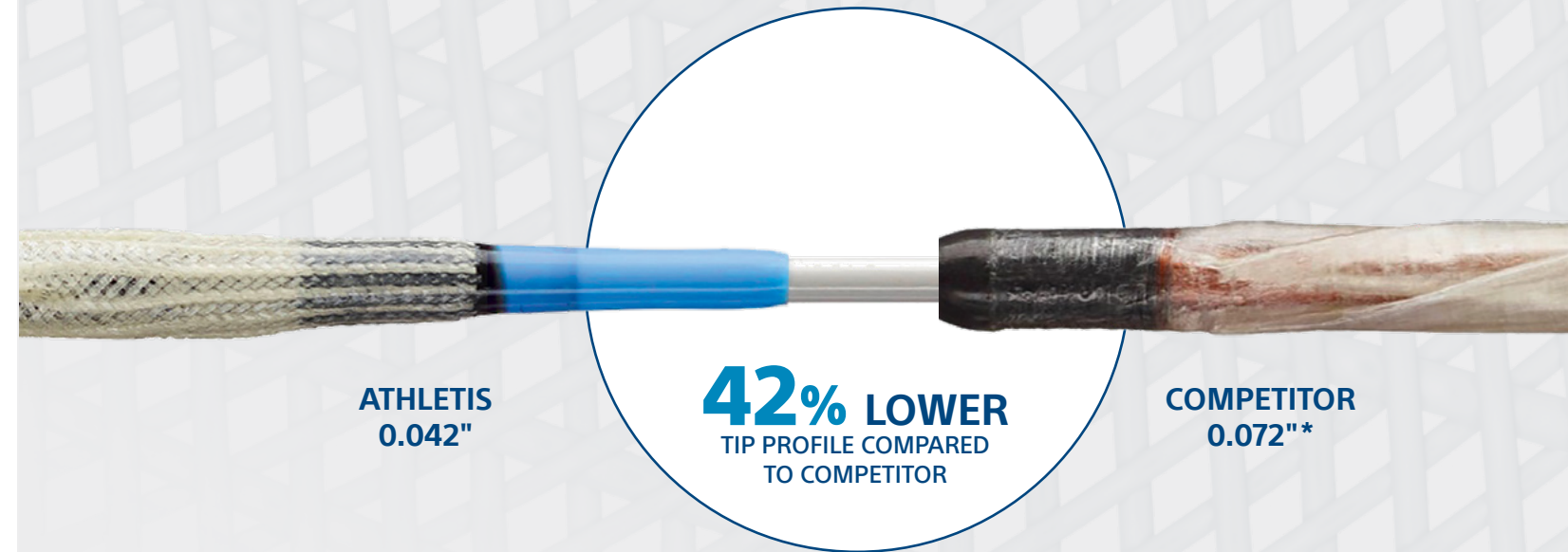

Athletis Ultra-High-Pressure Balloon Ordering Information

UPN GTIN			Balloon Length (mm)					
			20	30	40	60	80	100
6F	7F	8F						
Balloon Diameter (mm)	4	50	H74939347040250 08714729973539		H74939347040450 08714729973553	H74939347040650 08714729973577	H74939347040850 08714729973591	H74939347041050 08714729973614
		75	H74939347040270 08714729973546		H74939347040470 08714729973560	H74939347040670 08714729973584	H74939347040870 08714729973607	H74939347041070 08714729973621
	5	50	H74939347050250 08714729973638		H74939347050450 08714729973669	H74939347050650 08714729973683	H74939347050850 08714729973706	H74939347051050 08714729973720
		75	H74939347050270 08714729973645		H74939347050470 08714729973676	H74939347050670 08714729973690	H74939347050870 08714729973713	H74939347051070 08714729973737
		135			H74939347050410 08714729973652			
	6	50	H74939347060250 08714729973751		H74939347060450 08714729973782	H74939347060650 08714729973805	H74939347060850 08714729973829	H74939347061050 08714729973843
		75	H74939347060270 08714729973768		H74939347060470 08714729973799	H74939347060670 08714729973812	H74939347060870 08714729973836	H74939347061070 08714729973850
		135	H74939347060210 08714729973744		H74939347060410 08714729973775			
	7	50	H74939347070250 08714729973874		H74939347070450 08714729973904	H74939347070650 08714729973928	H74939347070850 08714729973942	H74939347071050 08714729973966
		75	H74939347070270 08714729973881		H74939347070470 08714729973911	H74939347070670 08714729973935	H74939347070870 08714729973959	H74939347071070 08714729973973
		135	H74939347070210 08714729973867		H74939347070410 08714729973898			
	8	50	H74939347080250 08714729973997	H74939347080350 08714729974017	H74939347080450 08714729974048	H74939347080650 08714729974062	H74939347080850 08714729974086	H74939347081050 08714729974109
		75	H74939347080270 08714729974000	H74939347080370 08714729974024	H74939347080470 08714729974055	H74939347080670 08714729974079	H74939347080870 08714729974093	H74939347081070 08714729974116
		135	H74939347080210 08714729973980		H74939347080410 08714729974031			
	9	50	H74939347090250 08714729974123		H74939347090450 08714729974154		H74939347090850 08714729974178	
		75	H74939347090270 08714729974130		H74939347090470 08714729974161		H74939347090870 08714729974185	
		135			H74939347090410 08714729974147			
	10	75	H74939347100270 08714729974192		H74939347100470 08714729974215		H74939347100870 08714729974222	
		135			H74939347100410 08714729974208			
	12	75	H74939347120270 08714729974246		H74939347120470 08714729974260	H74939347120670 08714729974284		
		135	H74939347120210 08714729974239		H74939347120410 08714729974253	H74939347120610 08714729974277		

INFLATION DEVICE ORDERING INFORMATION

Encore™ 40 Inflation Device is the recommended inflation device to be paired with Athletis Ultra-High-Pressure Balloon. To order, use **UPN M001394472540**.

*Average measurements taken by Boston Scientific. Data on file (TR 92200478). Bench test results may not necessarily be indicative of clinical performance.

ATHLETIS PTA BALLOON CATHETER

CAUTION: Federal law (USA) restricts this device to sale by or on the order of a physician. Rx only. Prior to use, please see the complete "Instructions for Use" for more information on Indications, Contraindications, Warnings, Precautions, Adverse Events, and Operator's Instructions. **INTENDED USE/ INDICATIONS FOR USE:** The Athletis PTA Balloon Dilatation Catheter is indicated for Percutaneous Transluminal Angioplasty in the peripheral vasculature including upper extremity, renal, iliac, and infrainguinal vessels and the treatment of obstructive lesions of native or synthetic arteriovenous dialysis fistulae. The Athletis PTA Balloon Dilatation Catheter is also indicated for post-dilatation of stents and stent grafts in the peripheral vasculature. **CONTRAINDICATIONS:** None known. **WARNINGS:** The Athletis PTA Balloon Dilatation Catheter is not for use in the coronary arteries. Any use for procedures other than those indicated in these instructions is not recommended. **PRECAUTIONS:** The Athletis™ PTA Balloon Dilatation Catheters are not intended for injection of contrast medium. To prevent over pressurization, use of a pressure monitoring device is recommended. Precautions to prevent or reduce clotting should be taken when any catheter is used: • Consider systemic anticoagulation. • Flush or rinse all products with sterile saline or a similar solution prior to use. **ADVERSE EVENTS:** The complications that may result from a balloon dilatation procedure include, but are not limited to: • Allergic reaction (device, contrast medium and medications or other) • Arrhythmia • Arteriovenous fistula • Cerebral Vascular Accident (CVA) / stroke / Transient Ischemic Attack (TIA) • Death • Embolization (air, device, plaque, thrombus or other) • Hematoma • Hemorrhage • Hypotension/hypertension • Ischemia • Myocardial Infarction/ischemia • Need for urgent intervention or surgery • Pain • Pulmonary Embolism • Renal Insufficiency or failure • Restenosis • Sepsis/infection • Thrombosis/thrombus • Vasospasm • Vessel injury (e.g. dissection, perforation, rupture, pseudoaneurysm, trauma, or other) • Vessel occlusion **92606536 A.1**

Athletis is a trademark of Boston Scientific Corporation or its affiliates. All other trademarks are property of their respective owners.

Boston Scientific
Advancing science for life™

Peripheral Interventions

300 Boston Scientific Way
Marlborough, MA 01752-1234
www.bostonscientific.com

To order product or for more information contact customer service at 1.888.272.1001.

© 2020 Boston Scientific Corporation or its affiliates. All rights reserved.

PI-810508-AB