

ZIPwire™

Nitinol Hydrophilic Guidewire

**Boston
Scientific**
Advancing science for life™

**The Fast Track for
Ureteral Access.**

Featuring
**Tip and Body Enhancements
to Stiff ZIPwire Guidewire**



ZIPwire™ Nitinol Hydrophilic Guidewire


The ZIPwire Guidewire is a kink-resistant nitinol wire with a lubricious hydrophilic coating designed for reliable ureteral access and reduced trauma. The fast track for ureteral access. The ZIPwire Guidewire is available in standard and stiff shaft.

Lubricious, floppy tip

facilitates reliable access to the ureter and reduced trauma during advancement through the urinary tract

Durable, covalently bonded hydrophilic coating designed to navigate ureteral strictures and facilitate tracking of urological instruments

Excellent torque control
when crossing tortuous anatomy

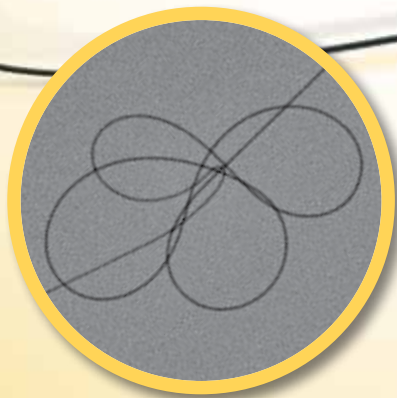

Highly radiopaque jacket

designed to provide enhanced visibility under fluoroscopy and enable accurate positioning



Torque vise may be used for added guidewire steerability and manipulation



Kink-resistant

designed for ease of advancement


The Stiffness Option

The Stiff ZIPwire Guidewire is designed to increase pushability, help straighten tortuous anatomy and facilitate advancement of devices.

Enhancements have been made to the tip and body of the current Stiff ZIPwire Guidewire to make it stiffer.

	Tip Enhancement	Body Enhancement
.025" Stiff ZIPwire Guidewire	✓	
.035" Stiff ZIPwire Guidewire	✓	✓
.038" Stiff ZIPwire Guidewire	✓	

According to bench testing, this table depicts the redesigned segments to be as rigid as other market devices. Our NEW Stiff ZIPwire Guidewires are designed to meet your stiffness needs.

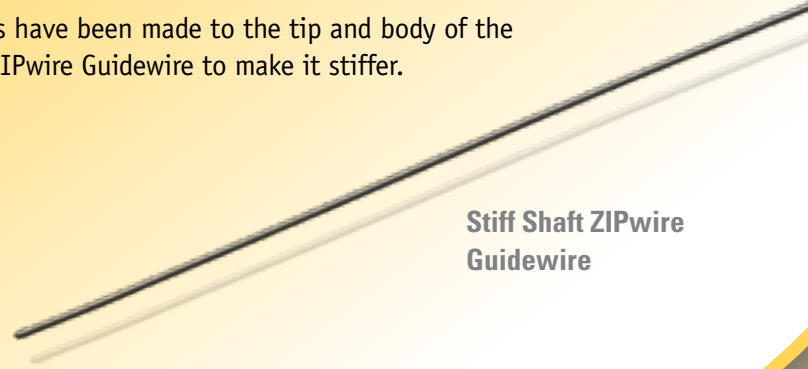


Stiff Shaft ZIPwire Guidewire

Multiple Tip Styles

Several tip styles are available, including the Bentson-type tip designed to provide the additional tip flexibility for negotiation of narrowed ureters and access beyond severe obstructions.

ZIPwire Guidewire Performance	
Reduced Trauma Soft, flexible tip designed for reduced trauma during advancement through the urinary tract	✓
Lubricious Hydrophilic coating designed for smooth advancement and easy tracking of devices	✓
Kink-Resistant Superelastic nitinol core designed to resist kinking	✓
Steerable Nitinol core and polyurethane jacket provide torqueability for excellent control and positioning	✓
Radiopaque Highly radiopaque jacket designed for enhanced fluoroscopic visibility	✓

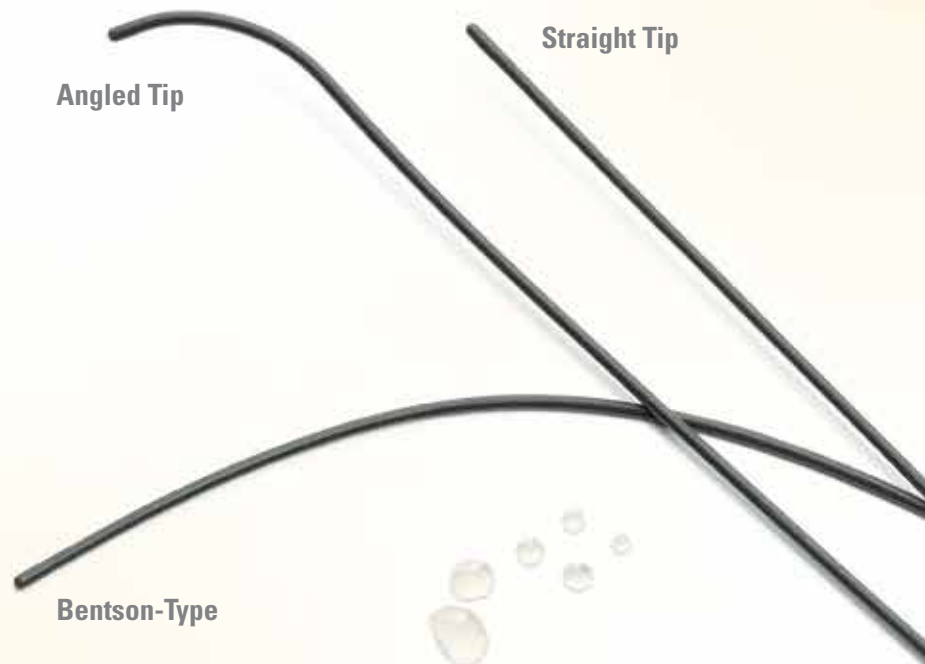

ORDER NUMBER	GTIN NUMBER	DIAMETER	LENGTH	TIP	SHAFT STIFFNESS	TAPER
Standard ZIPwire Guidewire						
M006630200B1	08714729802433	0.018 in	150cm	Straight	Standard	3cm
M006630201B1	08714729802457	0.025 in	150cm	Straight	Standard	3cm
M006630202B1	08714729802440	0.025 in	150cm	Angled	Standard	3cm
M006630212B1	08714729839149	0.032 in	150cm	Straight	Standard	3cm
M006630213B1	08714729839156	0.032 in	150cm	Angled	Standard	3cm
M006630205B1	08714729761761	0.035 in	150cm	Straight	Standard	3cm
M006630206B1	08714729761778	0.035 in	150cm	Angled	Standard	3cm
M006630208B1	08714729755333	0.038 in	150cm	Straight	Standard	3cm
M006630209B1	08714729761785	0.038 in	150cm	Angled	Standard	3cm
M006630210B1	08714729802556	0.038 in	260cm	Straight	Standard	3cm

Unit: Box 5

Stiff ZIPwire Guidewire						
M006630216B1	08714729839187	0.025 in	150cm	Straight	Stiff	3cm
M006630217B1	08714729839194	0.025 in	150cm	Angled	Stiff	3cm
M006630222B1	08714729755326	0.035 in	150cm	Straight	Stiff	3cm
M006630223B1	08714729761792	0.035 in	150cm	Angled	Stiff	3cm
M006630225B1	08714729761808	0.038 in	150cm	Straight	Stiff	3cm
M006630226B1	08714729802532	0.038 in	150cm	Angled	Stiff	3cm


Unit: Box 5

Bentson-Type ZIPwire Guidewire						
M006630214B1	08714729839163	0.035 in	150cm	Straight	Standard	Bentson 8cm
M006630203B1	08714729802464	0.035 in	150cm	Angled	Standard	Bentson 8cm
M006630215B1	08714729839170	0.038 in	150cm	Straight	Standard	Bentson 8cm
M006630207B1	08714729802501	0.038 in	150cm	Angled	Standard	Bentson 8cm

Unit: Box 5

Bentson-Type Stiff ZIPwire Guidewire						
M006630224B1	08714729802495	0.035 in	150cm	Straight	Stiff	Bentson 5cm
M006630221B1	08714729802488	0.035 in	150cm	Angled	Stiff	Bentson 5cm
M006630227B1	08714729802549	0.038 in	150cm	Straight	Stiff	Bentson 5cm
M006630228B1	08714729802525	0.038 in	150cm	Angled	Stiff	Bentson 5cm

Unit: Box 5

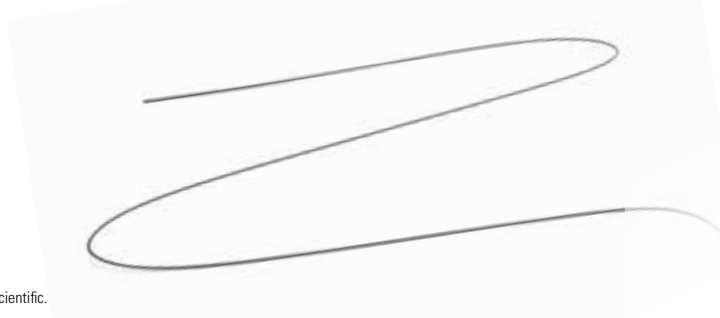

Boston Scientific
Advancing science for life™

Boston Scientific Corporation
300 Boston Scientific Way
Marlborough, MA 01752-1234
www.bostonscientific.com

Ordering Information
1.888.272.1001

© 2016 Boston Scientific Corporation
or its affiliates. All rights reserved.

URO-135002-AC NOV 2016

All images are the property of Boston Scientific.

All trademarks are the property of their respective owners.

Products shown for INFORMATION purposes only and may not be approved or for sale in certain countries. Please check availability with your local sales representative or customer service.

Caution: U.S. Federal law restricts this device to sale by or on the order of a physician.

CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labelling supplied with each device. Information for use only in countries with applicable health authority registrations. Material not intended for use in France.