

AccuStick™ II

Introducer Systems With Radiopaque Marker

Boston
Scientific

For the introduction of
.038" guidewires in
non-vascular procedures

RADIOPAQUE MARKER

RO marker positioned at tip of sheath enhances visualization for more precise placement

LARGE SHEATH END HOLE INNER DIAMETER

Facilitates placement of .038" working guidewire alongside .018" guidewire

CONVENIENT, INTEGRATED SYSTEM


- ◆ 21-Gauge diagnostic needle reduces tissue trauma
- ◆ Option of .018" kink-resistant nitinol guidewire is designed to provide strength during tough procedures
- ◆ Coaxial sheath/dilator assembly with locking stiffening cannula is designed for over-the-wire placement
- ◆ .038" heavy duty guidewire in a choice of "J" tip or straight tip provides excellent performance as a working wire

SHEATH/DILATOR ASSEMBLY PROXIMAL TAPER

Eases introduction and placement

DILATOR REFERENCE MARK


Facilitates alignment of sheath and dilator tips for more accurate wire placement

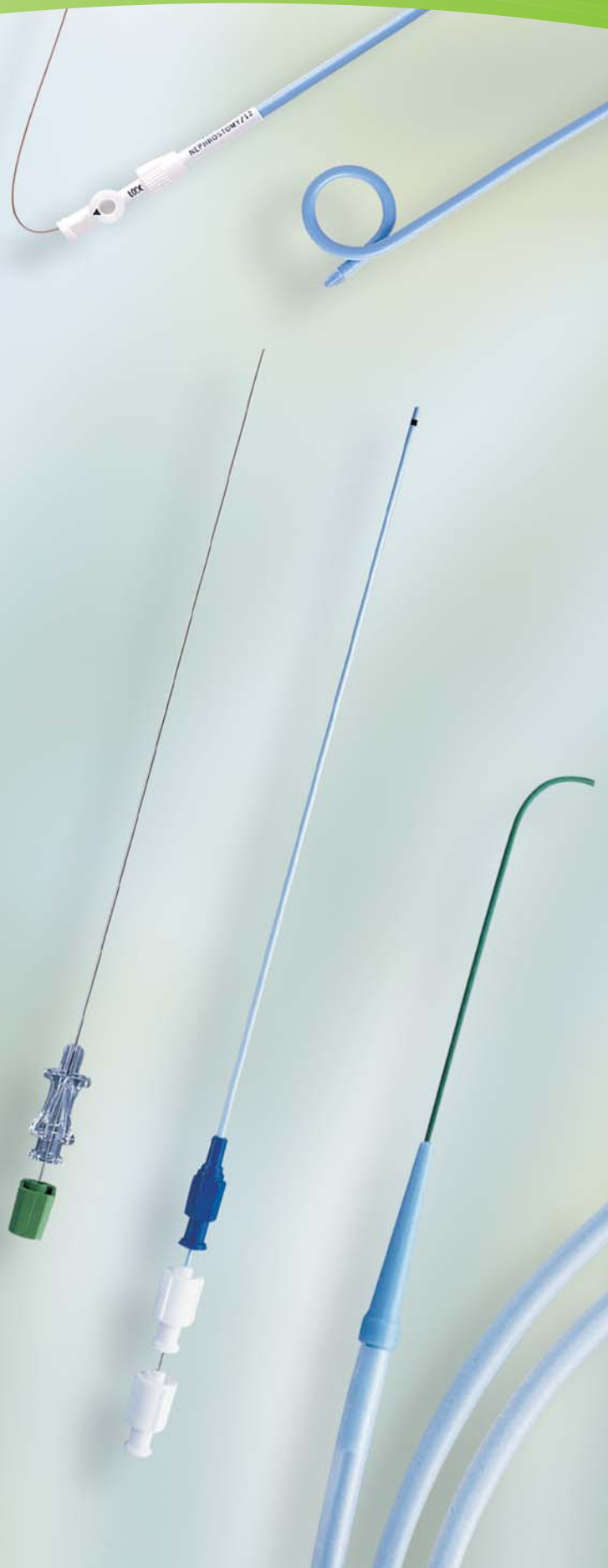


AccuStick™ II

Introducer Systems With Radiopaque Marker


The .018" guidewire is advanced fully through the needle.

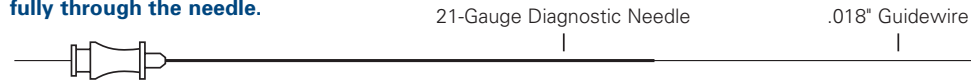

After removing the needle, the locked sheath and dilator/cannula assembly is advanced over the .018" guidewire.

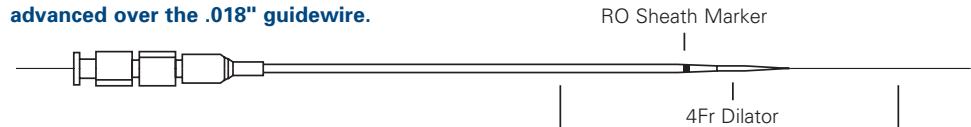

The dilator/cannula assembly is removed and the .038" guidewire is advanced through the 6Fr sheath.

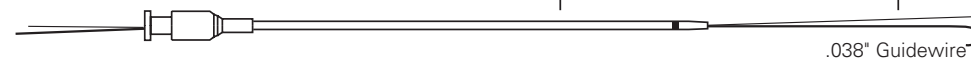

ORDERING INFORMATION

UPN	ORDER NUMBER	DESCRIPTION	QUANTITY
M001207020*	20-702	AccuStick II Introducer System with - 1 Extra Stiff, Stainless Steel, 60cm, .018" Guidewire with 7.5cm Floppy Platinum Tip - 1 Heavy Duty, Stainless Steel, 145cm, .038" Guidewire with 3mm "J" Tip	1
M001207030*	20-703	AccuStick II Introducer System with - 1 Extra Stiff, Stainless Steel, 60cm, .018" Guidewire with 7.5cm Floppy Platinum Tip - 1 Heavy Duty, Stainless Steel, 145cm, .038" Guidewire with Straight Tip	1
M001207100*	20-710	AccuStick II Introducer System with - 1 Nitinol, 60cm, .018" Guidewire with 6cm Floppy Platinum Tip - 1 Heavy Duty, Stainless Steel, 145cm, .038" Guidewire with 3mm "J" Tip	1
M001207040*	20-704	AccuStick II Introducer System with - 1 Extra Stiff, Stainless Steel, 60cm, .018" Guidewire with 7.5cm Floppy Platinum Tip	1
M001207050*	20-705	AccuStick II Introducer System without Guidewires	1
M001201081	20-108	Stainless Steel, 60cm, .018" Guidewire with 7.5cm Floppy Platinum Tip	5
M001207071	20-707	Stainless Steel, 60cm, .018" Guidewire with 11.5cm Floppy Platinum Tip	5
M001207111	20-711	Nitinol, 60cm, .018" Guidewire with 6cm Floppy Platinum Tip	5
M001206041	20-604	21-Gauge Diagnostic Needle with Stylet	10

*Each system includes: One AccuStick II Introducer (Coaxial 6Fr Sheath with RO Marker/4Fr Dilator, 20cm) with Locking Stiffening Cannula; One 21-Gauge, 15cm Diagnostic Needle with Stylet

INDICATIONS: The AccuStick™ II Introducer System with Radiopaque Marker facilitates introduction and placement of a guidewire.

CONTRAINDICATIONS: None known

PRECAUTIONS: Advancement and/or withdrawal of .018"/.038" (0.46mm/0.97mm) wires should be smooth and without resistance. **Do not use excessive force.** If resistance is felt, carefully remove wire(s) and system as a unit. Device should be used for non-vascular procedures.

CAUTION: Federal (USA) law restricts this device to sale by or on the order of a physician. Refer to package insert provided with the product for complete Instructions for Use, Contraindications, Potential Adverse Effects, Warnings and Precautions prior to using this product.

Boston Scientific

Delivering what's next.™

Boston Scientific Corporation
Oncology
One Boston Scientific Place
Natick, MA 01760-1537
www.bostonscientific.com

To order product or for more information, contact Customer Service at **1.800.225.3228**

© 2005 by Boston Scientific Corporation or its affiliates. All rights reserved.

eONCPS170 / 5M / 02/05