


AGENT™

Drug-Coated Balloon

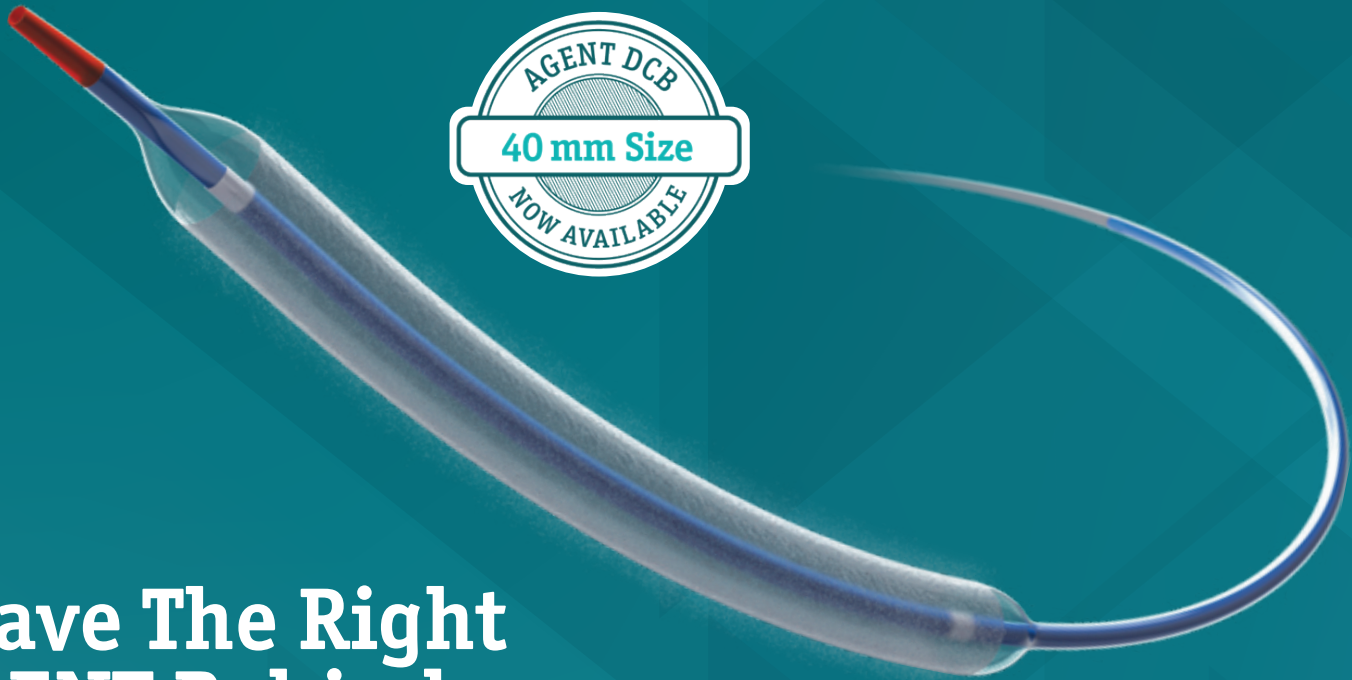

Leave The Right AGENT Behind

AGENT makes Reliable Drug Transfer More Deliverable.

Compliance Chart

Pressure atm (kPa)	Balloon Size (mm)						
	2.00	2.25	2.50	2.75	3.00	3.50	4.00
3.0 (304)	1.86	2.06	2.28	2.53	2.76	3.19	3.66
4.0 (405)	1.93	2.14	2.37	2.61	2.85	3.30	3.80
5.0 (507)	1.99	2.20	2.44	2.68	2.93	3.39	3.88
6.0 (608)	2.03	2.26	2.50	2.75	3.00	3.46	3.96
7.0 (709)	2.07	2.31	2.55	2.81	3.06	3.52	4.04
8.0 (811)	2.10	2.34	2.59	2.85	3.11	3.57	4.14
9.0 (912)	2.13	2.38	2.62	2.88	3.15	3.61	4.00
10.0 (1013)	2.15	2.40	2.65	2.91	3.18	3.64	4.18
11.0 (1115)	2.18	2.42	2.67	2.94	3.21	3.68	4.22
12.0 (1216)	2.19	2.44	2.69	2.96	3.23	3.72	4.25
13.0 (1317)	2.21	2.46	2.72	2.99	3.26	-	-
14.0 (1419)	2.23	2.48	2.74	3.02	3.28	-	-


Nominal Pressure Rated Burst Pressure. Do Not Exceed.

*Drug coating designed to deliver therapeutic drug dose in association with balloon angioplasty treatment.

† Bench testing as compared to Apex™ and Maverick™ balloons. Data on file at Boston Scientific.

AGENT Drug-Coated Balloon vis distributed by Boston Scientific Corporation. The legal manufacturer is Hemoteq AG, Adenauerstrasse 15, 52146 Würselen Germany.

AGENT™ Drug-Coated Balloon

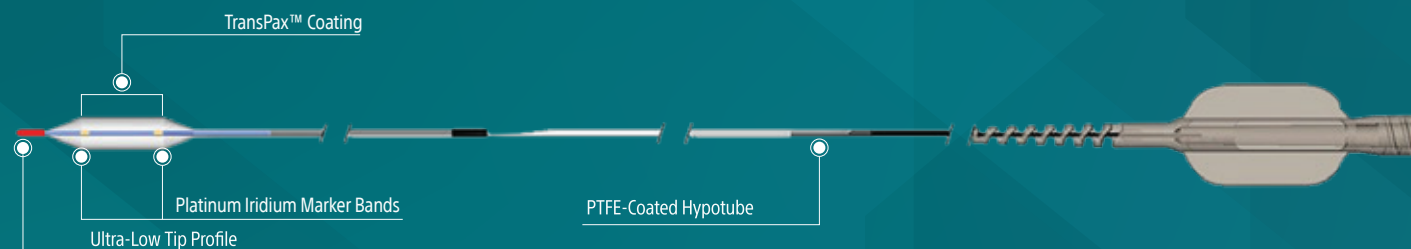

Product Details

Indications for Use	The AGENT DCB Balloon Catheter is indicated for Percutaneous Transluminal Coronary Angioplasty (PTCA) in the coronary arteries to treat in-stent restenosis (ISR) and de novo / small vessel disease.
Drug Coating	The TransPax drug coating is a proprietary formulation of Paclitaxel and a citrate ester excipient (acetyl tributyl citrate - ATBC). The AGENT balloon catheter is designed to prevent hyperplasia of smooth muscle cells at the treated lesion site.
Drug Dose	2.0 µg per mm² of the balloon surface
Recommended Inflation Time	At least 30 seconds
Antiplatelet Regimen	Recommended for 3 months after treatment
Available Balloon Lengths	In mm: 12, 15, 20, 30, 40
Available Balloon Diameters	In mm: 2.00, 2.25, 2.50, 2.75, 3.00, 3.50, 4.00
Lesion Entry Profile	0.43 mm or 0.017"
Delivery System Effective Length	144 cm
Delivery System Port	Designed for 0.36 mm or 0.014" guidewire
Balloon Material	OptiLEAP™
Shelf Life	24 months for Use
Sterilization	Ethylene oxide (EO)
Single/Multiple Use	Single Use
Guide Catheter Inner Diameter	1.42 mm or 0.056"
Guide Catheter Compatibility	5 F or 1.67 mm
Catheter Shaft Outer Diameter	1.67 mm or 0.0656"
Marker Band Material	Platinum-Iridium
Marker Band Placement	The outside edges of the marker bands indicate the balloon's shoulders and coated region

Ordering Information

Balloon Length (mm)					
ø (mm)	12mm	15mm	20mm	30mm	40mm
2.00	H749 3922220121 0	H749 3922220151 0	H749 3922220201 0	H749 3922220301 0	H749 3922220401 0
2.25	H749 3922222121 0	H749 3922222151 0	H749 3922222201 0	H749 3922222301 0	
2.50	H749 3922225121 0	H749 3922225151 0	H749 3922225201 0	H749 3922225301 0	H749 3922225401 0
2.75	H749 3922227121 0	H749 3922227151 0	H749 3922227201 0	H749 3922227301 0	
3.00	H749 3922230121 0	H749 3922230151 0	H749 3922230201 0	H749 3922230301 0	H749 3922230401 0
3.50	H749 3922235121 0	H749 3922235151 0	H749 3922235201 0	H749 3922235301 0	H749 3922235401 0
4.00	H749 3922240121 0	H749 3922240151 0	H749 3922240201 0	H749 3922240301 0	H749 3922240401 0

**Boston
Scientific**

Advancing science for life™

Boston Scientific India Private Limited
Corporate Office: 3rd Floor, Bestech
Business Towers, Sector-48, Sohna Road,
Gurugram - 122018, India
Phone - +91-124-6260500
Registered Address: Unit No 304, 3rd Floor,
Chiranjiv Tower, 43 Nehru Place,
New Delhi 110019
For customer queries:
customercare.india@bsci.com
www.bostonscientific.com/en-IN

The information presented in this document by Boston Scientific is intended for healthcare professionals in India. [The information presented in this document is intended to provide balanced, scientific, and evidence based answers to unsolicited medical questions.] The information should not be, however, treated as comprehensive and does not intend to provide diagnosis, treatment or any medical or health advice. For full information about the product, please refer to the appropriate product labeling. Responsibility for patient care resides with the healthcare professional on the basis of his or her professional licence, experience and knowledge of the patient. Healthcare professionals must rely on their judgment when deciding which treatments and procedures to use with patients.

All cited trademarks are the property of their respective owners.
©2026 Boston Scientific Corporation or its affiliates. All rights reserved.

INDH-IC-2498401-May 2026