

Dreamtome™ RX

Sphincterotome

**Boston
Scientific**
Advancing science for life™



RX Biliary System™

We're in!™

Dreamtome™ RX

Sphincterotome

High Performance... Simplified

Sphincterotome pre-loaded with Dreamwire™ Guidewire

- Facilitates a smooth transition from over-the-wire to short-wire ERCP technique without removing devices from the endoscope
- 260cm and 450cm guidewire lengths available to meet physician preference for short-wire or traditional over-the-wire technique
- Striped movement markers designed to indicate wire positioning during wire advancement and device exchange

RX Biliary System™

We're in!™



The RX Biliary System helps put the physician in control of the guidewire. Using the RX Locking Device, the physician can secure the guidewire to maintain access while exchanging devices.


The RX Biliary System and its compatible devices offer a broad portfolio of products for pancreatico-biliary stone and stricture management.

Guidewire Design and Tip Technology

Dream Access:

- 10cm hydrophilic Dream Tip™ End is designed for soft, smooth, access and excellent tactile feel
- Endoglide™ Coating intended to promote smooth tracking and enhance tactile sensation
- Guidewire shaft designed for efficient device exchanges
- Forms an Alpha Loop to facilitate effective access into the targeted duct

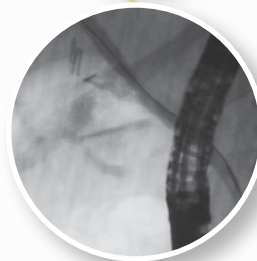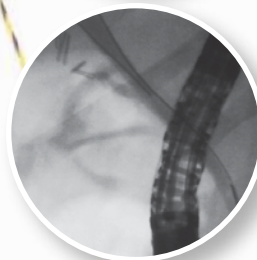


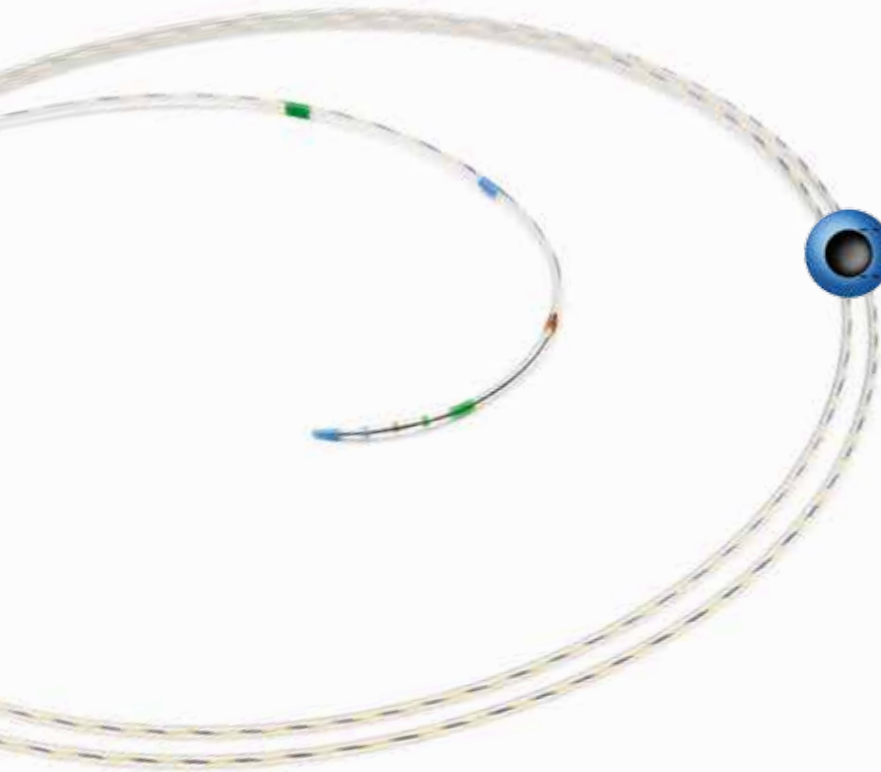

Merging Lumen Technology™

Proximal lumens merge into a single distal lumen on RX 44 Sphincterotome

- Developed to facilitate contrast injection while maintaining guidewire access through open lumen technology

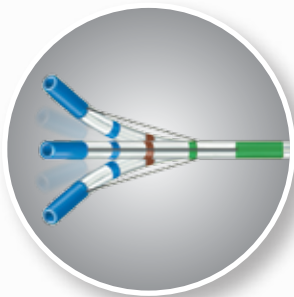

Flexible and Responsive Design


Direct Wire Technology™ is designed for flexibility and responsiveness to aid in precise sphincterotomy

- Tip orientation and cut-wire positioning capabilities developed to facilitate cannulation
- Ergonomic, rotating handle engineered for excellent control
- Tapered tip sizes to accommodate varying patient anatomy

RX LOCKING DEVICE[†]

Guidewire locks in place to maintain access


- Available with short 260cm guidewire or a standard 450cm guidewire to meet physician technique preferences
- Guidewire may be locked and unlocked at any time during the procedure
- Facilitates device exchange, contrast injection and biliary therapy

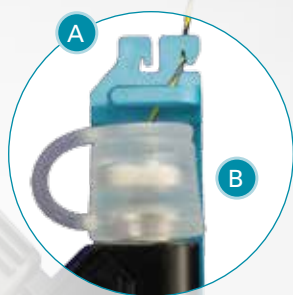


- A** Locking device has notches for fixation of one or two wires
- B** Biopsy port cap designed to reduce air and bile leakage
- C** Strap secures it below the biopsy port

[†]RX Locking Device sold separately.

DREAM ACCESS

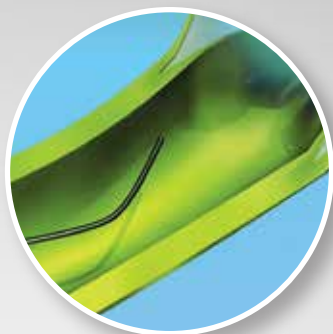

Dreamwire™ Guidewire gains access into targeted duct

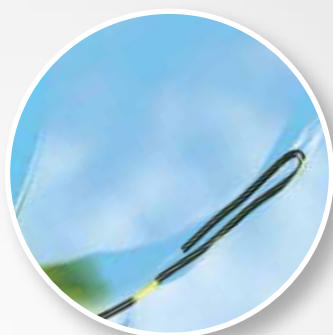


Dreamwire Guidewire encounters an obstacle (stone, stricture) and forms an Alpha Loop



Dreamwire Guidewire successfully passes through obstacle, enabling the next step in the ERCP

Dreamtome™ RX

Sphincterotome

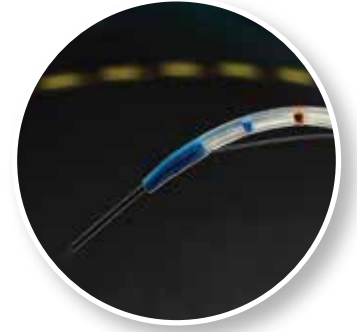

Dreamtome RX Sphincterotome

Order Number	GTIN	Description	Cut Wire (mm)	Tip Length (mm)	Preloaded Guidewire (in / cm)	Guidewire Tip Style	Packaging
M00584000	08714729778226	Dreamtome RX 49 Sphincterotome	20	5	.035 / 260 Dreamwire Guidewire	Straight	Each
M00584010	08714729778233	Dreamtome RX 49 Sphincterotome	30	5	.035 / 260 Dreamwire Guidewire	Straight	Each
M00584020	08714729778240	Dreamtome RX 49 Sphincterotome	20	5	.035 / 450 Dreamwire Guidewire	Straight	Each
M00584030	08714729778257	Dreamtome RX 49 Sphincterotome	30	5	.035 / 450 Dreamwire Guidewire	Straight	Each
M00584040	08714729778264	Dreamtome RX 44 Sphincterotome	20	5	.035 / 260 Dreamwire Guidewire	Straight	Each
M00584050	08714729778271	Dreamtome RX 44 Sphincterotome	30	5	.035 / 260 Dreamwire Guidewire	Straight	Each
M00584060	08714729778288	Dreamtome RX 44 Sphincterotome	20	5	.035 / 450 Dreamwire Guidewire	Straight	Each
M00584070	08714729778295	Dreamtome RX 44 Sphincterotome	30	5	.035 / 450 Dreamwire Guidewire	Straight	Each

Preloaded with straight-tip guidewires only. Active Cords sold separately.

All photographs taken by Boston Scientific

CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labelling supplied with each device. Information for use only in countries with applicable health authority registrations. This material not intended for use in France. RX Only.

All trademarks are the property of their respective owners.

Boston Scientific
Advancing science for life™

Boston Scientific Corporation
300 Boston Scientific Way
Marlborough, MA 01752-1234
www.bostonscientific.com

Ordering Information
1.888.272.1001

©2019 Boston Scientific Corporation or its affiliates. All rights reserved.

ENDO-377607-AB