

# Dreamwire™

High Performance Guidewire

**Boston  
Scientific**

Advancing science for life™




*RX Biliary System™*

**We're in!™**

## High Performance Guidewire


### Dream Access - *Facilitates deep cannulation of targeted ducts*

- 10cm hydrophilic Dream Tip™ End is designed for soft, smooth, access and excellent tactile feel
- Choice of angled or straight Dream Tip End
- Endoglide™ Coating promotes smooth tracking and enhances tactile sensation
- Forms an Alpha Loop to facilitate effective access into the targeted duct

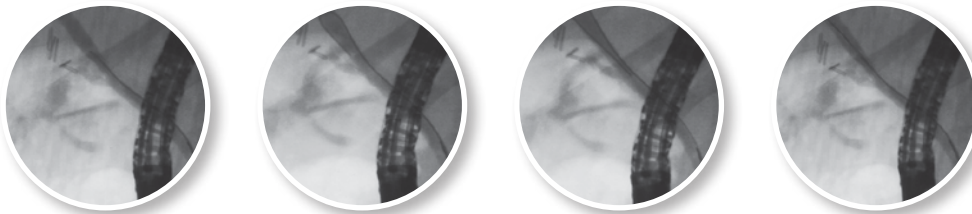

### Dream Exchange - *Designed to improve procedural efficiency and easy RX Biliary device*

- 260cm and 450cm guidewire lengths available to meet physician preference for short-wire or traditional over-the-wire technique
- Kink-resistant Nitinol standard and stiff shaft, engineered for wire control and manipulation capabilities
- Striped movement markers designed to indicate wire positioning during wire advancement and device exchange
- Side loading Torque Device designed for effective guidewire tip manipulation (Angled-tip only)
- Guidewire may be left in place during sphincterotomy

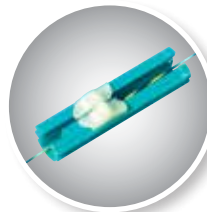

### RX Biliary System™

# We're in!™

The RX Biliary System helps put the physician in control of the guidewire. Using the RX Locking Device, the physician can secure the guidewire to maintain access while exchanging devices.

The RX Biliary System and its compatible devices offer a broad portfolio of products for pancreatico-biliary stone and stricture management.

## Dreamwire High Performance Guidewire

Order Number	GTIN	Description	Distal Tip O.D. (in)	Working Length (cm)	Shaft Style	Dream Tip Style	Tip
M00556101	08714729777816	Dreamwire Guidewire ST	0.035	260	Standard	Straight	Round
M00556111	08714729777847	Dreamwire Guidewire ANG	0.035	260	Standard	Angled	Round
M00556121	08714729777830	Dreamwire Guidewire ST Stiff	0.035	260	Stiff	Straight	Round
M00556131	08714729777823	Dreamwire Guidewire ANG Stiff	0.035	260	Stiff	Angled	Round
M00556141	08714729777885	Dreamwire Guidewire ST	0.035	450	Standard	Straight	Round
M00556151	08714729777878	Dreamwire Guidewire ANG	0.035	450	Standard	Angled	Round
M00556161	08714729777861	Dreamwire Guidewire ST Stiff	0.035	450	Stiff	Straight	Round
M00556171	08714729777854	Dreamwire Guidewire ANG Stiff	0.035	450	Stiff	Angled	Round

ST = Straight, ANG = Angled, Packaging: Box 2

All trademarks are the property of their respective owners.

Indications, Contraindications, Warnings and Instructions for Use can be found in the product labeling supplied with each device.

**Caution:** Federal (USA) law restricts this device to sale by or on the order of a physician.

## Boston Scientific

Advancing science for life™

Boston Scientific Corporation  
300 Boston Scientific Way  
Marlborough, MA 01752-1234  
www.bostonscientific.com

**Ordering Information**  
**1.888.272.1001**

©2016 Boston Scientific Corporation  
or its affiliates. All rights reserved.

ENDO-376410-AA May 2016