

Standard Percutaneous Endoscopic Gastrostomy (PEG) Kit

Standard PEG Kit

Order Number	Description	Units
M00568381	20Fr (6.67mm) Standard PEG Kit, Pull	Box 2
M00568401	20Fr (6.67mm) Standard PEG Kit, Push	Box 2
M00568391	24Fr (8.00mm) Standard PEG Kit, Pull	Box 2
M00568411	24Fr (8.00mm) Standard PEG Kit, Push	Box 2

PEG Accessories

Order Number	Description	Units
M00580641	20Fr (6.67mm) Y-Ports	Box 10
M00580651	24Fr (8.00mm) Y-Ports	Box 10
M00580731	20Fr (6.67mm) Round Bolsters	Box 10
M00580781	20Fr (6.67mm) and 24Fr (8.00mm) C-Clamps	Box 10

Through-The-PEG (TTP) Jejunal Feeding tube

EndoVive TTP Jejunal Feeding Tube

Order Number	Description	Tip Style	Units
M00566320	8.5Fr (2.83mm) TTP Jejunal Feeding TUBE	Bent Tip	Box 1
M00566330	8.5Fr (2.83mm) TTP Jejunal Feeding TUBE	Pigtail Tip	Box 1
M00566340	12Fr (4.00mm) TTP Jejunal Feeding TUBE	Bent Tip	Box 1
M00566350	12Fr (4.00mm) TTP Jejunal Feeding TUBE	Pigtail Tip	Box 1

Recommended Anchoring Device

Order Number	Description	Working Length (cm)	Minimum Working Channel (mm)	Clip Opening (mm)	Units
M00522610	Resolution™ Clip	235	2.8	11	Each
M00522611	Resolution™ Clip	235	2.8	11	Box 10
M00522612	Resolution™ Clip	235	2.8	11	Box 20

Gastrostomy Tubes

Standard Profile Balloon Replacement Feeding Tubes

Order Number	Description	Units
M00582040	12Fr (4.00mm) Replacement G-Tube, Straight	Box 1
M00582050	14Fr (4.67mm) Replacement G-Tube, Straight	Box 1
M00582060	16Fr (5.33mm) Replacement G-Tube, Straight	Box 1
M00582070	18Fr (6.00mm) Replacement G-Tube, Straight	Box 1
M00582080	20Fr (6.67mm) Replacement G-Tube, Straight	Box 1
M00582090	22Fr (7.33mm) Replacement G-Tube, Straight	Box 1
M00582100	24Fr (8.00mm) Replacement G-Tube, Straight	Box 1
M00582110	12Fr (4.00mm) Replacement G-Tube, Right-Angle	Box 1
M00582120	14Fr (4.67mm) Replacement G-Tube, Right-Angle	Box 1
M00582130	16Fr (5.33mm) Replacement G-Tube, Right-Angle	Box 1
M00582140	18Fr (6.00mm) Replacement G-Tube, Right-Angle	Box 1
M00582150	20Fr (6.67mm) Replacement G-Tube, Right-Angle	Box 1
M00582161	22Fr (7.33mm) Replacement G-Tube, Right-Angle	Box 1
M00582171	24Fr (8.00mm) Replacement G-Tube, Right-Angle	Box 1

Low Profile Button Replacement Feeding Tubes

Low Profile Button Replacements

Order Number	Description	Length (cm)	Color Code	Units
M00562760	18Fr (6.00mm) Low Profile Button Replacement	1.2	Yellow	Box 1
M00562770	18Fr (6.00mm) Low Profile Button Replacement	1.7	Yellow	Box 1
M00562780	18Fr (6.00mm) Low Profile Button Replacement	2.4	Yellow	Box 1
M00562790	18Fr (6.00mm) Low Profile Button Replacement	3.4	Yellow	Box 1
M00562800	18Fr (6.00mm) Low Profile Button Replacement	4.4	Yellow	Box 1
M00562820	24Fr (8.00mm) Low Profile Button Replacement	1.7	Blue	Box 1
M00562830	24Fr (8.00mm) Low Profile Button Replacement	2.4	Blue	Box 1
M00562840	24Fr (8.00mm) Low Profile Button Replacement	3.4	Blue	Box 1
M00562850	24Fr (8.00mm) Low Profile Button Replacement	4.4	Blue	Box 1
M00568170	28Fr (9.33mm) Low Profile Button Replacement	1.5	Green	Box 1
M00568180	28Fr (9.33mm) Low Profile Button Replacement	2.8	Green	Box 1
M00568190	28Fr (9.33mm) Low Profile Button Replacement	4.3	Green	Box 1

Low Profile Button Accessories

Order Number	Description	Length (cm)	Color Code	Units
M00580401	18Fr (6.00mm) Decompression Tube	1.2	Yellow	Box 10
M00580411	18Fr (6.00mm) Decompression Tube	1.7	Yellow	Box 10
M00580421	18Fr (6.00mm) Decompression Tube	2.4	Yellow	Box 10
M00580431	18Fr (6.00mm) Decompression Tube	3.4	Yellow	Box 10
M00580441	18Fr (6.00mm) Decompression Tube	4.4	Yellow	Box 10
M00580461	24Fr (8.00mm) Decompression Tube	1.7	Blue	Box 10
M00580471	24Fr (8.00mm) Decompression Tube	2.4	Blue	Box 10
M00580481	24Fr (8.00mm) Decompression Tube	3.4	Blue	Box 10
M00580491	24Fr (8.00mm) Decompression Tube	4.4	Blue	Box 10
M00580171	18Fr (6.00mm) Bolus Straight Adapter		Yellow	Box 10
M00580181	18Fr (6.00mm) Continuous Right-Angle Adapter		Yellow	Box 10
M00580231	24Fr (8.00mm) Bolus Straight Adapter		Blue	Box 10
M00580241	24Fr (8.00mm) Continuous Right-Angle Adapter		Blue	Box 10
M00580271	28Fr (9.33mm) Bolus Straight Feeding Adapter		Green	Box 10
M00580281	28Fr (9.33mm) Continuous Right-Angle Feeding Adapter		Green	Box 10

EndoVive® Enteral Access Devices

Boston
Scientific



**Boston
Scientific**
Advancing science for life™

www.bostonscientific.eu

All cited trademarks are the property of their respective owners. CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labelling supplied with each device. Information for the use only in countries with applicable health authority product registrations.

©2016 Boston Scientific Corporation or its affiliates. All rights reserved.

ENDO-301501-AA February 2016.

Product Features

High-Grade Silicone Dome-Bolstered PEG Tube

- Designed for performance and consistent flow
- Facilitates traction removal
- Universal 1.9 mm O.D. x 240 cm retrieval snare

Universal 1.9 mm O.D. x 240 cm Retrieval Snare

- Traverses scopes with a 2.0 mm or larger working channel
- Universal-sized snare, intended to eliminate the cost of a specialty snare

Multiple Sizes and Accessories


- Multiple French sizes for broad patient coverage
- Separate feeding ports
- C-Clamp designed to help reduce risk of tube reflux and leaking

Double-Tier Tray Packaging

- Kit organisation promotes procedural efficiency

External Bolsters

- Two external bolster options: accommodate physician and patient needs: Round; T-Bar

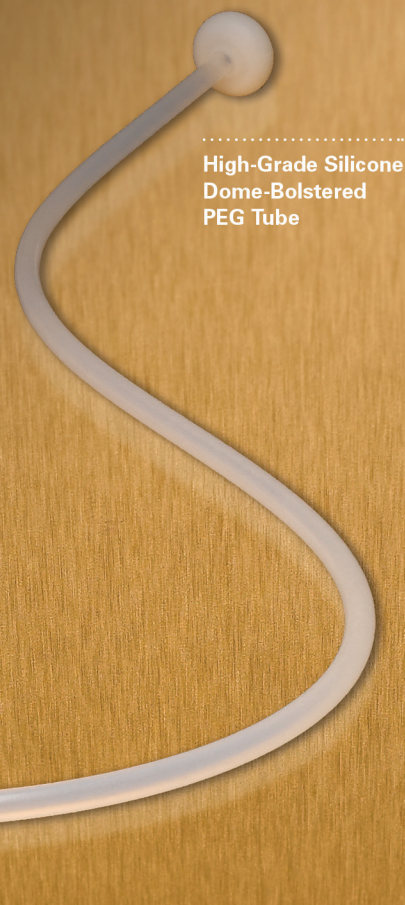

High-Grade Silicone Dome-Bolstered PEG Tube

Through-The-PEG (TTP) Jejunal Feeding Tube

Product Features

Innovative Device Design

- Designed for placement through an existing Standard PEG
 - 68 cm length, 8.5Fr TTP J-Tube compatible with 20Fr PEG
 - 68 cm length, 12Fr TTP J-Tube compatible with 24Fr PEG

- Available in Pigtail Tip and Bent Tip designs, each with a suture loop intended to facilitate distal tip anchoring past the Ligament of Treitz

- Compatible with both Pull and Push (over-the-wire) Placement Techniques

Flexima Material

- Radiopaque material facilitates fluoroscopic visualisation

Temp Tip


- Temporary tapered Tip designed to provide smooth tracking into the jejunum to help ensure proper placement
- Temp Tip is designed to dissolve at body temperature permitting utilisation of the full distal tip lumen diameter which may help reduce premature tube clogging

Hollow Stiffening Stylet


- Provides added stiffness designed to improve pushability and permit greater physician control during placement

Recommended Anchoring Device


- Resolution™ Clip

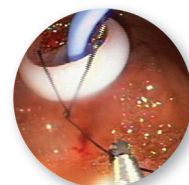

Resolution™ Clip grasping

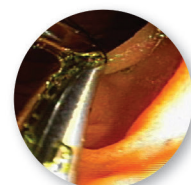

Resolution™ Clip anchoring suture and jejunal tube to gastric wall

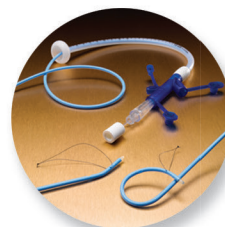

TTP Jejunal Feeding Tube



Resolution™ Clip


Product Features

“Pancake” Shaped Balloon

- Designed to offer the deflation characteristics of a balloon for insertion and removal that may be more comfortable for the patient
- Consistent 6 ml balloon inflation across all sizes

Standard Profile Gastrostomy with Optional Right-Angled External Bolster

- Designed to minimise the potential tube kinking when it's taped down
- Designed to minimise premature pull out due to shape of the tube
- Engineered to minimise the likelihood that the Y-Port and the feeding tube adapter will come apart and spill the contents of the nutrition bag
- Separate medication and feeding ports
- Clear shaft for internal lumen visualisation

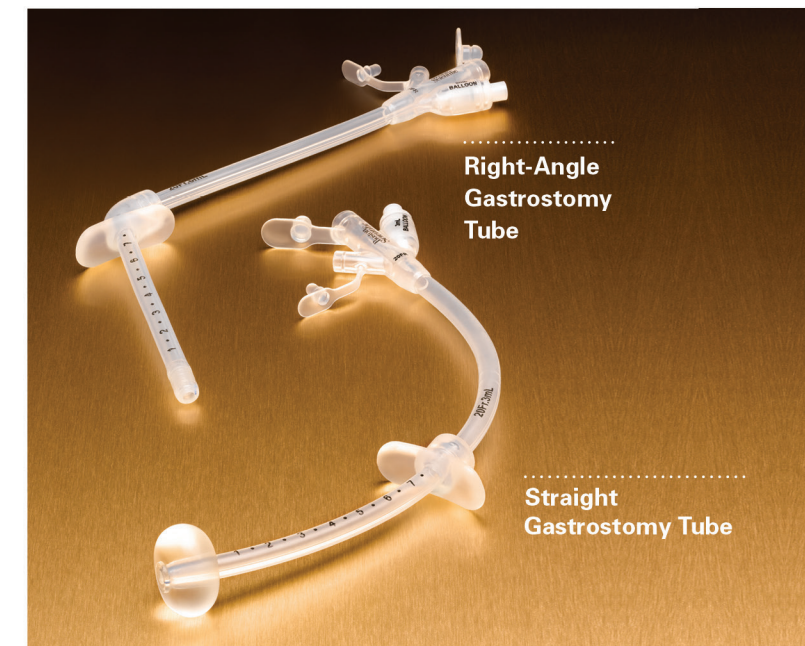

Right-Angle Gastrostomy Tube

Straight Gastrostomy Tube

KIT CONTENTS

- (1) Standard Balloon Replacement with Straight or Right-Angle External Bolster
- (2) 4"x 4" Gauze Pads
- (1) 6 ml Syringe
- (1) Lubricating Jelly

Low Profile Button Replacements


Product Features

High-Grade Silicone Button-Bolstered Low Profile Replacement

- Radiopaque silicone designed to enhance placement confirmation
- Inserted through an existing stoma tract
- Low profile design helps reduce risk of accidental pull-out
- Product flexibility at skin level designed to enhance patient comfort and cosmetic appeal
- Designed to reduce gastric backflow
- Kit includes accessories designed to allow for gastric decompression
- Color-coded French sizes with corresponding feeding accessories and multiple lengths for broad patient coverage
- Specially designed round spacers to accommodate “in-between” sizes and short-term weight gain

KIT CONTENTS

- (1) Button Kit
- (2) Percutaneous Stoma Measuring Device (PSMD)
- (1) Spacer Disk Variety Pack
- (1) Right Angle Decompression Tube
- (1) Patient Care Kit

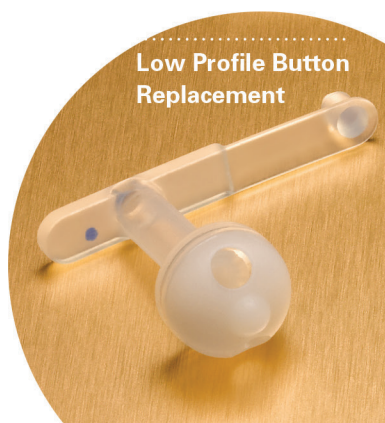

Low Profile Button Replacement