

Jagwire™ Revolution


High Performance Guidewire

GET WIRED FOR SUCCESS

To advance and manage pancreaticobiliary diseases efficiently and effectively you need a suite of access and cannulation solutions, supported by clinical evidence and designed to deliver the best patient care possible.

Successful cannulation is a critical first step for all ERCP procedures and having the right guidewire can make a real difference, from the most simple to the most complex cases.

**Introducing the Jagwire Revolution
High Performance Guidewire...**

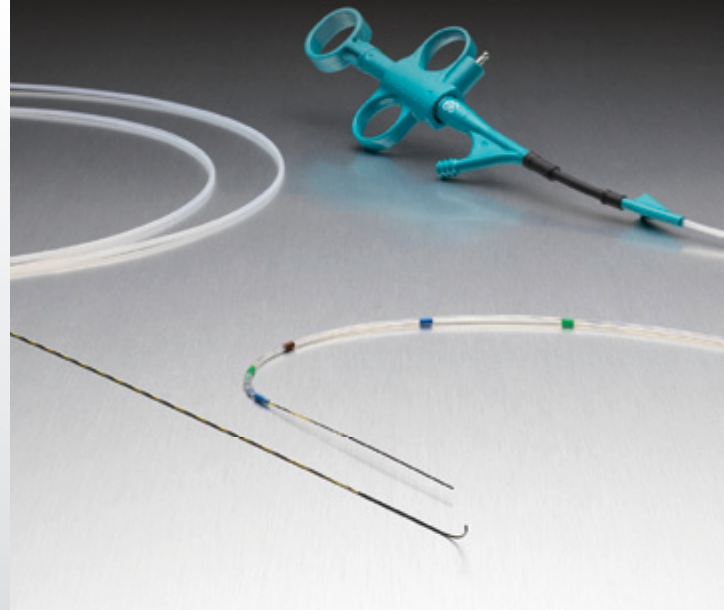

Similar Stiffness and Pushability of our .035" in a .025" Guidewire

The Jagwire Revolution High Performance Guidewire is a versatile .025" guidewire engineered to have the stiffness and pushability that you would expect in our Jagwire .035" guidewire combined with both access and exchange characteristics to help improve completion rates of initial ERCPs.

Efficient Access and Navigation

- Tight alpha-loop formation designed for knuckling through challenging ducts and strictures
- Atraumatic 5 cm hydrophilic tapered tip
- Super elastic Nitinol core designed to enhance performance in tortuous anatomy

Navigation in Tortuous Anatomy


- Hydrophilic coating designed to aid passability through devices and difficult strictures
- Straight and angled tips for physician preference and difficult anatomy
- Engineered to facilitate precise torque ability of the angled tip

Improved Visualization Design


- Tungsten filled tip for improved fluoroscopic visualization
- Striped guidewire movement markers

Control and Durability

- RX Biliary System compatible design helps enable clinician to maintain control of the guidewire throughout the procedure and engage the locking mechanism
- Innovative and durable coating design

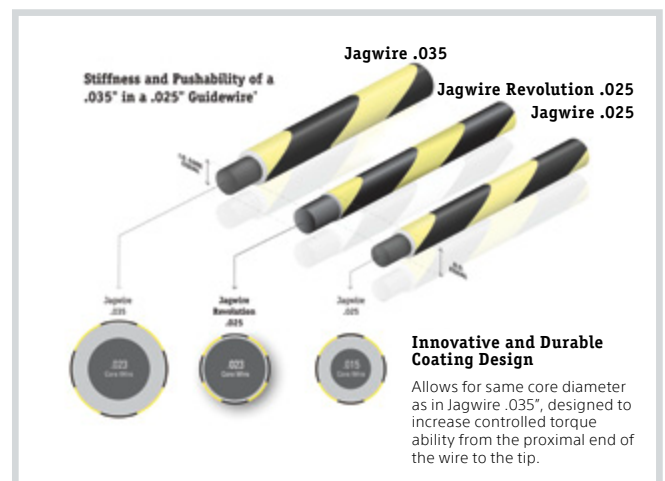

*Data on File PDM E0445

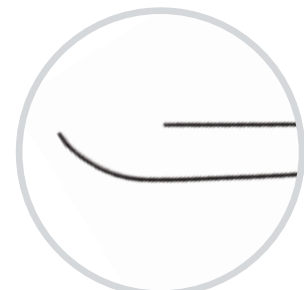

Indications for Use

The Jagwire Revolution High Performance Guidewire is indicated for use in selective cannulation of the biliary ducts including the common bile duct, pancreatic duct, cystic duct, and right and left hepatic ducts. The endoscopic guidewire is designed to be used during endoscopic pancreaticobiliary procedures for catheter introduction, exchanges of catheters, cannulas, and sphincterotomes, and to aid in the placement of diagnostic and therapeutic devices.

Ordering Information: Jagwire Revolution High Performance Guidewire

Order Number	Description	Distal Tip O.D. (in)	Working Length (cm)	Tip Style	Tip	Packaging (e.g., Box 2)
M00557001	Jagwire Revolution High Performance Guidewire- ST	.025	260	Straight	Round	Box 2
M00557011	Jagwire Revolution High Performance Guidewire- ANG	.025	260	Angled	Round	Box 2
M00557021	Jagwire Revolution High Performance Guidewire- ST	.025	450	Straight	Round	Box 2
M00557031	Jagwire Revolution High Performance Guidewire- ANG	.025	450	Angled	Round	Box 2
M00557041	Jagwire Revolution High Performance Guidewire- ST	.025	500	Straight	Round	Box 2

Ordering Information: Jagtome™ Revolution RX Cannulating Sphincterotome (Autotome™ Pre-loaded with Jagwire Revolution High Performance Guidewire)

Order Number	Description	Cut Wire (mm)	Tip Length (mm)	Distal Tip O.D. (F)	Pre-loaded Guidewire (in/cm)	Guidewire Tip Style	Packaging (e.g., each)
M00584240	Jagtome Revolution RX 39 Sphincterotome	20	5	3.9	.025/260 Jagwire Revolution Guidewire	Straight	Each
M00584250	Jagtome Revolution RX 39 Sphincterotome	20	5	3.9	.025/450 Jagwire Revolution Guidewire	Straight	Each
M00584260	Jagtome Revolution RX 39 Sphincterotome	30	5	3.9	.025/260 Jagwire Revolution Guidewire	Straight	Each
M00584270	Jagtome Revolution RX 39 Sphincterotome	30	5	3.9	.025/450 Jagwire Revolution Guidewire	Straight	Each

Ordering Information: TRUEtome™ Jagwire™ Revolution Cannulating Sphincterotome (TRUEtome™ Pre-loaded with Jagwire Revolution High Performance Guidewire)

Order Number	Description	Cut Wire (mm)	Tip Length (mm)	Distal Tip O.D. (F)	Pre-loaded Guidewire (in/cm)	Guidewire Tip Style	Packaging (e.g., each)
M00583260	TRUEtome Jagwire Revolution 39 Sphincterotome	20	5	3.9	.025/450 Jagwire Revolution Guidewire	Straight	Each
M00583280	TRUEtome Jagwire Revolution 39 Sphincterotome	30	5	3.9	.025/450 Jagwire Revolution Guidewire	Straight	Each

Boston Scientific Contact Information

Name: _____ Title: _____

Phone Number: _____ E-mail: _____

The testing was performed by or on behalf of Boston Scientific. Data on file. All photographs taken by Boston Scientific. Products shown for INFORMATION purposes only and may not be approved or for sale in certain countries. Please check availability with your local sales representative or customer service.

All cited trademarks are the property of their respective owners. CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labeling supplied with each device. Information for the use only in countries with applicable health authority product registrations.

ENDO-595706-AB Printed in Germany by medicalvision.

Boston Scientific
Advancing science for life™

www.bostonscientific.eu

© 2019 Boston Scientific Corporation or its affiliates. All rights reserved.
DINEND2409EB