

# Hurricane™ RX

## Biliary Balloon Dilatation Catheter

### HURRICANE RX SINGLE-USE BILIARY BALLOON DILATATION CATHETER\*

Order Number	Balloon O.D. (mm)	Balloon Length (cm)	Catheter O.D. (Fr) / (mm)	Catheter Length (cm)	Maximum Inflation Pressure (ATM) / (kPa)	Minimum Channel Size (mm)	Recommended Guidewire (in) / (mm)
M00545890	4	2	5.8 / 1.93	180	12 / 1216	3.2	.035" / 0.89
M00545900	4	4	5.8 / 1.93	180	12 / 1216	3.2	.035" / 0.89
M00545910	6	2	5.8 / 1.93	180	11 / 1114	3.2	.035" / 0.89
M00545920	6	4	5.8 / 1.93	180	11 / 1114	3.2	.035" / 0.89
M00545930	8	2	5.8 / 1.93	180	11 / 1114	3.2	.035" / 0.89
M00545940	8	4	5.8 / 1.93	180	11 / 1114	3.2	.035" / 0.89
M00545950	10	2	5.8 / 1.93	180	8 / 810	3.2	.035" / 0.89
M00545960	10	4	5.8 / 1.93	180	8 / 810	3.2	.035" / 0.89

\*Recommended Guidewire: .035" (0.89mm) Hydra Jagwire™ Guidewire.

### ENCORE™ 26 INFLATION DEVICE

M00566670 Encore 26 Inflation Device

**European Distribution Centre – The Netherlands**  
T: +31 45 54 67 700 F: +31 45 54 67 800  
**European Headquarters – Paris**  
T: +33 1 57 66 80 00 F: +33 1 57 66 84 99  
**Argentina (Freephone)**  
T: 0800 555 2678 F: +5411 4896 8517  
**Australia/New Zealand**  
T: +61 2 8063 8100 F: +61 2 9330 1404  
**Austria**  
T: +43 1 60 810 F: +43 1 60 810 60  
**Belgium (Freephone)**  
T: 0800 94 494 F: 0800 93 343  
**Brazil**  
T: +55 11 5502 8500 F: +55 11 5103 2212  
**Canada**  
T: +1 888 359 9691 F: +1 888 575 7396  
**Chile**  
T: +562 445 4904 F: +562 445 4915  
**China – Beijing**  
T: +86 10 8525 1588 F: +86 10 8525 1566  
**China – Guangzhou**  
T: +86 20 8767 9791 F: +86 20 8767 9789  
**China – Shanghai**  
T: +86 21 6141 5959 F: +86 21 6141 5900

**Colombia**  
T: +57 1 629 5045 F: +57 1 612 4761  
**Czech Republic**  
T: +420 296 331 901 F: +420 296 331 935  
**Denmark (Freephone)**  
T: 80 30 80 02 F: 80 30 80 05  
**Eire (Freephone)**  
T: 1800 882 969 F: 1800 882 968  
**Finland**  
T: +358 20 762 8882 F: +358 20 762 8883  
**France**  
T: +33 1 39 30 49 00 F: +33 1 39 30 49 01  
**Germany (Freephone)**  
T: 08000 723300 F: 08000 723319  
**Greece**  
T: +30 210 9542 300 F: +30 210 9542 310  
**Hong Kong**  
T: +852 2960 7100 F: +852 2563 5276  
**Hungary**  
T: +36 1 456 30 40 F: +36 1 456 30 41  
**India – Bangalore**  
T: +91 80 2212 4928/9 F: +91 80 2207 5153  
**India – Chennai**  
T: +91 44 2220 1879 F: +91 44 2220 0716  
**India – Delhi**  
T: +91 11 4243 2222 F: +91 11 2610 0808

**India – Mumbai**  
T: +91 22 4030 9165 F: +91 22 4040 9199  
**Italy**  
T: +39 010 60 60 1 F: +39 010 60 60 200  
**Korea**  
T: +82 2 3476 2121 F: +82 2 3476 1776  
**Mexico**  
T: +52 55 5687 63 90 F: +52 55 5687 62 28  
**Middle East/Gulf/North Africa**  
T: +961 1 805 410 F: +961 1 805 445  
**The Netherlands**  
T: +31 30 602 55 44 F: +31 30 602 55 05  
**Norway (Freephone)**  
T: 800 104 04 F: 800 101 90  
**Poland**  
T: +48 22 435 14 14 F: +48 22 435 14 10  
**Portugal**  
T: +351 1 381 25 40 F: +351 21 381 25 58  
**South Africa**  
T: +27 11 840 8600 F: +27 11 463 6077  
**South East Asia – Malaysia**  
T: +60 3 2263 3813 F: +60 3 2264 3813  
**South East Asia – Philippines**  
T: +63 2 687 6994 F: +63 2 687 3047  
**South East Asia – Singapore**  
T: +65 6418 8888 F: +65 6418 8889

**South East Asia – Thailand**  
T: +66 2 6380 100 F: +66 2 6380 400  
**South East Europe**  
T: +30 210 9542 300 F: +30 210 9542 310  
**Spain – Barcelona**  
T: +34 93 444 72 00 F: +34 93 405 90 45  
**Spain – Madrid**  
T: +34 901 11 12 15 F: +34 91 319 50 03  
**Sweden**  
T: +46 42 25 69 00 F: +46 42 25 69 69  
**Switzerland**  
T: +41 32 626 57 00 F: +41 32 626 57 01  
**Taiwan**  
T: +886 2 2747 7278 F: +886 2 2747 7270  
**Turkey – Istanbul**  
T: +90 216 464 36 66 F: +90 216 464 36 67  
**UK (Freephone)**  
T: 0844 800 4512 F: 0844 800 4513  
**Uruguay**  
T: +59 82 900 6212 F: +59 82 900 6212  
**Venezuela**  
T: +58 212 959 6275 F: +58 212 959 5328

All cited trademarks are the property of their respective owners. CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labelling supplied with each device. Information for the use only in countries with applicable health authority product registrations.  
PSST 3118 Printed in the UK by Gosling.

**Boston Scientific**  
Delivering what's next.™

www.bostonscientific.com  
www.bostonscientific-international.com

© 2009 Boston Scientific Corporation  
or its affiliates. All rights reserved.  
DINEND2125ED

# Hurricane™ RX

## Biliary Balloon Dilatation Catheter

Boston  
Scientific



RX Biliary System™  
**A difference you can feel.™**

# RX Biliary System™

A difference you can feel.™

The endoscope shows you the way, but it's your sense of touch that gets you where you need to go. The devices that make up the RX Biliary System have been developed to provide access and control with a sensation of device responsiveness that truly is... *a difference you can feel.*


Over the years, Boston Scientific has earned the confidence and allegiance of customers around the world who have adopted the RX Biliary System for everyday use. Our legacy of success in developing high performance RX biliary devices has enabled us to change the landscape of ERCP procedures. By innovating and refining products, we've put control and access right at your fingertips.

## Hurricane™ RX Biliary Balloon Dilatation Catheter



### ACHIEVES PRESSURE UP TO 12ATM (1216kPa)

- ◆ Long 4cm balloon length is designed to offer stability during stricture dilation
- ◆ Semi-compliant material engineered to adapt to biliary duct anatomy during balloon inflation
- ◆ Balloon is recommended for use with the Encore™ 26 Inflation Device which has been developed for rapid deflation to potentially reduce patient discomfort

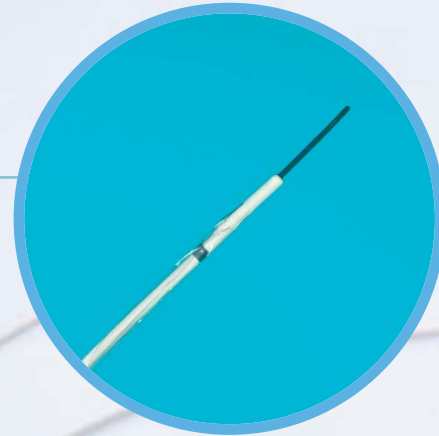

### EXCELLENT STRICTURE PASSAGE

*Low-profile, 5.8Fr (1.93mm) balloon shaft is designed to reduce resistance and increase pushability when crossing strictures*

- ◆ 4Fr (1.33mm) tip is designed for less traumatic advancement across the stricture
- ◆ Stainless steel stiffening stylet inserted in the proximal portion of the catheter shaft to augment pushability
- ◆ RX catheter design allows guidewire locking to secure access

### DESIGNED FOR ACCURATE POSITIONING

*Hurricane RX Balloon offers multiple options for visualisation*

- ◆ Two radiopaque markers aid in balloon positioning
- ◆ Contrast may be injected with guidewire in place



Dilation of a stricture with the Hurricane Dilatation Balloon

### RANGE OF BALLOON SIZES

*Options of balloon diameters and lengths facilitate selection of appropriately sized balloon for patient anatomy and stricture size*

- ◆ 2cm and 4cm lengths available
- ◆ 4mm, 6mm, 8mm and 10mm balloon diameters available