



# Single-use Broncho Videoscope with Full HD Visualization Endoscopic Image Processor

**Boston  
Scientific**  
Advancing science for life™







# Scope Features



## Three distinct sizes of Bronchoscope

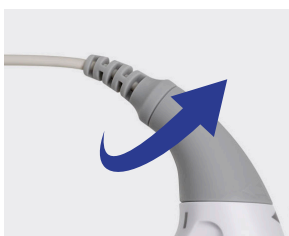

120° rotation of rear end of operating handle

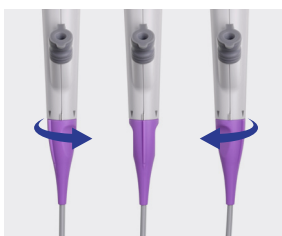

The  $\pm 60^\circ$  shaft rotation supports insertion of scope during procedures



The direction of the suction button is adjustable, which can prevent tangled pipes



Two-way bending 220° Up/Down



Passability of accessory devices



A wide and clear field of vision with 120° field of view and 0.5-120mm dept of field

The Bronchus Sampler is designed as an accessory to the Single-use Broncho Videoscope. It is intended to enable the aspiration and collection of fluid samples from the bronchi or alveoli areas of the lungs with up to 60ml of collection per sample container.



No Reprocessing



No Cross Infection



No Maintenance

The product is a sterile single-use device, and the device must be used in a manner consistent with IFU.



## Full HD visualization endoscopic image processor



\*Some features may require a USB drive for storage which is not included with device



# Specifications



## Single-Use Broncho Videoscope

Model	SBV-1A-B	SBV-1B-B	SBV-1C-B
Field of View	120°	120°	120°
Depth of View Field	0.5-120mm	0.5-120mm	0.5-120mm
Outer Diameter of Insertion Section (OD)	2.8mm	4.2mm	5.6mm
Working Channel (ID)	1.2mm	2.0mm	2.8mm
Working Length	600mm	600mm	600mm
Total Length	920mm	920mm	920mm
Angulation Range	220° Up/Down	220° Up/Down	220° Up/Down



## Full HD Visualization Endoscopic Image Processor HDVS-S300D

Size	332mm (W) x 240mm (H) x 50mm (D)
Screen Size (mm)	13.3" TFT-LCD
Resolution	1920 x 1080p
Output	HDMI x 1, SDI x 1
IPX	IPX0
Weight	2.5Kg
USB Interface	Type-A x 1, Type-C x 1
Storage	External USB Storage
Battery Operation	4H

## Product Ordering Codes

Product Code	Description	Unit
SBV-1A-B	Single-use Broncho Videoscope, Distal End OD: 2.7mm, Insertion Tube OD: 2.8mm	Ea
SBV-1B-B	Single-use Broncho Videoscope, Distal End OD: 4.0mm, Insertion Tube OD: 4.2mm	Ea
SBV-1C-B	Single-use Broncho Videoscope, Distal End OD: 5.6mm, Insertion Tube OD: 5.6mm	Ea
BS-6060	Bronchus Sampler	Box 5
SC-60	Sample Container	Box 10
BRK-S01	Bracket	Ea
HDVS-S300D	Full HD Visualization Endoscopic Image Processor HDVS-S300D	Ea

Name of registrant/manufacture: Scivita Medical Technology (Suzhou) Co., Ltd.

CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings, and instructions for use can be found in the product labelling supplied with each device for both Boston Scientific and Scivita Medical products, or at [www.IFU-BSCI.com](http://www.IFU-BSCI.com) for applicable Boston Scientific's products. Products shown for INFORMATION purposes only and may not be approved or for sale in certain countries. This material not intended for use in France. All trademarks are property of their respective owners

Exclusively distributed by

**Boston Scientific**

Advancing science for life™

©2025 Boston Scientific Corporation or its affiliates. All rights reserved. ENDO-2218301-AC