

Trapezoid™ RX

Wire-guided Retrieval Basket

TRAPEZOID RX WIRE-GUIDED RETRIEVAL BASKET

Order Number	Description	Open Basket Diameter (cm)	Required Working Channel (mm)	Recommended Guidewire* (in) / (mm)
M00510860	Trapezoid RX	1.5	3.2	0.035 / 0.89
M00510870	Trapezoid RX	2.0	3.2	0.035 / 0.89
M00510880	Trapezoid RX	2.5	3.2	0.035 / 0.89
M00510890	Trapezoid RX	3.0	3.2	0.035 / 0.89

*Recommended guidewire: 0.035" (0.89mm) Hydra Jagwire™ Guidewire

ALLIANCE™ II INFLATION SYSTEM

Order Number	Description
M00550620	Inflation Handle (each)

European Distribution Centre – The Netherlands
T: +31 45 54 67 700 F: +31 45 54 67 800

European Headquarters – Paris
T: +33 1 57 66 80 00 F: +33 1 57 66 84 99

Argentina (Freephone)
T: 0800 555 2678 F: +5411 4896 8517

Australia/New Zealand
T: +61 2 8063 8100 F: +61 2 9330 1404

Austria
T: +43 1 60 810 F: +43 1 60 810 60

Belgium (Freephone)
T: 0800 94 494 F: 0800 93 343

Brazil
T: +55 11 5502 8500 F: +55 11 5103 2212

Canada
T: +1 888 359 9691 F: +1 888 575 7396

Chile
T: +562 445 4904 F: +562 445 4915

China – Beijing
T: +86 10 8525 1588 F: +86 10 8525 1566

China – Guangzhou
T: +86 20 8767 9791 F: +86 20 8767 9789

China – Shanghai
T: +86 21 6141 5959 F: +86 21 6141 5900

Colombia
T: +57 1 629 5045 F: +57 1 612 4761

Czech Republic
T: +420 296 331 901 F: +420 296 331 935

Denmark (Freephone)
T: 80 30 80 02 F: 80 30 80 05

Eire (Freephone)
T: 1800 882 969 F: 1800 882 968

Finland
T: +358 20 762 8882 F: +358 20 762 8883

France
T: +33 1 39 30 49 00 F: +33 1 39 30 49 01

Germany (Freephone)
T: 0800 723300 F: 0800 723319

Greece
T: +30 210 9542 300 F: +30 210 9542 310

Hong Kong
T: +852 2960 7100 F: +852 2563 5276

Hungary
T: +36 1 456 30 40 F: +36 1 456 30 41

India – Bangalore
T: +91 80 2212 4928/9 F: +91 80 2207 5153

India – Chennai
T: +91 44 2220 1879 F: +91 44 2220 0716

India – Delhi
T: +91 11 4243 2222 F: +91 11 2610 0808

India – Mumbai
T: +91 22 4030 9165 F: +91 22 4040 9199

Italy
T: +39 010 60 60 1 F: +39 010 60 60 200

Korea
T: +82 2 3476 2121 F: +82 2 3476 1776

Mexico
T: +52 55 5687 63 90 F: +52 55 5687 62 28

Middle East/Gulf/North Africa
T: +961 1 805 410 F: +961 1 805 445

The Netherlands
T: +31 30 602 55 44 F: +31 30 602 55 05

Norway (Freephone)
T: 800 104 04 F: 800 101 90

Poland
T: +48 22 435 14 14 F: +48 22 435 14 10

Portugal
T: +351 1 381 25 40 F: +351 21 381 25 58

South Africa
T: +27 11 640 9800 F: +27 11 463 6077

South East Asia – Malaysia
T: +60 3 2263 3813 F: +60 3 2284 3813

South East Asia – Philippines
T: +63 2 687 6994 F: +63 2 687 3047

South East Asia – Singapore
T: +65 6418 8888 F: +65 6418 8899

South East Asia – Thailand
T: +66 2 6380 100 F: +66 2 6380 400

South East Europe
T: +30 210 9542 300 F: +30 210 9542 310

Spain – Barcelona
T: +34 93 444 72 00 F: +34 93 405 90 45

Spain – Madrid
T: +34 901 11 12 15 F: +34 91 319 50 03

Sweden
T: +46 42 25 69 00 F: +46 42 25 69 69

Switzerland
T: +41 32 626 57 00 F: +41 32 626 57 01

Taiwan
T: +886 2 2747 7278 F: +886 2 2747 7270

Turkey – Istanbul
T: +90 216 464 36 66 F: +90 216 464 36 67

UK (Freephone)
T: 0844 800 4512 F: 0844 800 4513

Uruguay
T: +59 82 900 6212 F: +59 82 900 6212

Venezuela
T: +58 212 959 6275 F: +58 212 959 5328

Boston Scientific

Delivering what's next.™

www.bostonscientific.com
www.bostonscientific-international.com

All cited trademarks are the property of their respective owners. CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labelling supplied with each device. Information for the use only in countries with applicable health authority product registrations.
PSST IC 1828 Printed in the UK by Gosling.

© 2008 Boston Scientific Corporation
or its affiliates. All rights reserved.
DINEND2153EB

Trapezoid™ RX

Wire-guided Retrieval Basket

Boston Scientific



RX Biliary System™
A difference you can feel.™

RX Biliary System™

A difference you can feel.™

The endoscope shows you the way, but it's your sense of touch that gets you where you need to go. The devices that make up the RX Biliary System have been developed to provide access and control with a sensation of device responsiveness that truly is... *a difference you can feel.*

Over the years, Boston Scientific has earned the confidence and allegiance of customers around the world who have adopted the RX Biliary System for everyday use. Our legacy of success in developing high performance RX biliary devices has enabled us to change the landscape of ERCP procedures. By innovating and refining products, we've put control and access right at your fingertips.

DESIGNED FOR OPTIMAL PERFORMANCE

Flexi-Stiff™ Catheter Construction and metal alloy basket provide a distinct difference in stone retrieval

- ◆ Flexible and sturdy, Flexi-Stiff Catheter design developed to bend with biliary anatomy and maintain excellent pushability and tactile sensation
- ◆ Metal alloy basket material is engineered to maintain basket symmetry, even after difficult stone manipulation




MECHANICAL STONE CRUSHING OPTION

Alliance™ II Inflation Device engineered to crush stones with scope in place


- ◆ The inflation device can be used for mechanical advantage to crush large, difficult stones trapped within the basket

Trapezoid™ RX

Wire-guided Retrieval Basket

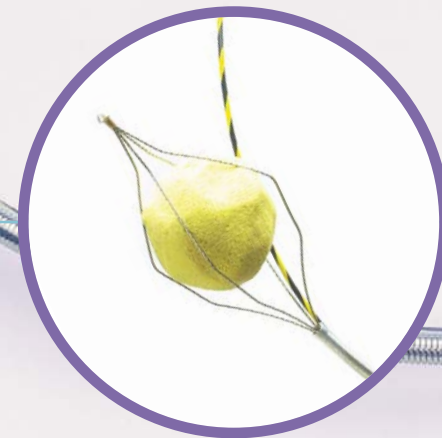

- ◆ Tip disengages to facilitate basket removal if stone does not break, reducing risk of basket impaction

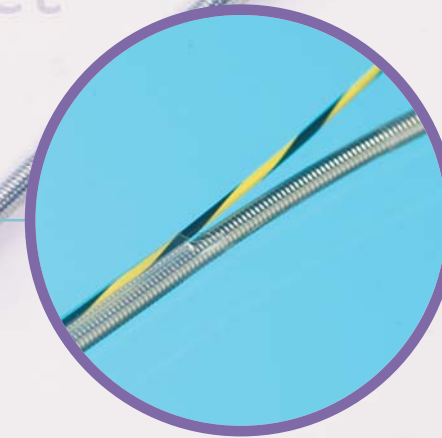

MAINTAINS GUIDEWIRE ACCESS

Distal side-port guidewire channel facilitates duct access

- ◆ Side-port designed to lead basket past stones and help retrieve stones
- ◆ Distal side-port channel accepts 0.035" (0.89mm) short 260cm guidewire or 450cm traditional guidewire to accommodate physician technique preference
- ◆ Injection with wire in place accommodates cholangiography without device exchange
- ◆ Wire locking device facilitates multiple sweeps of the duct during stone retrieval by maintaining guidewire position
- ◆ Wire can be locked and unlocked at any time during the procedure