


# Hydra Jagwire™

Guidewire

Boston  
Scientific

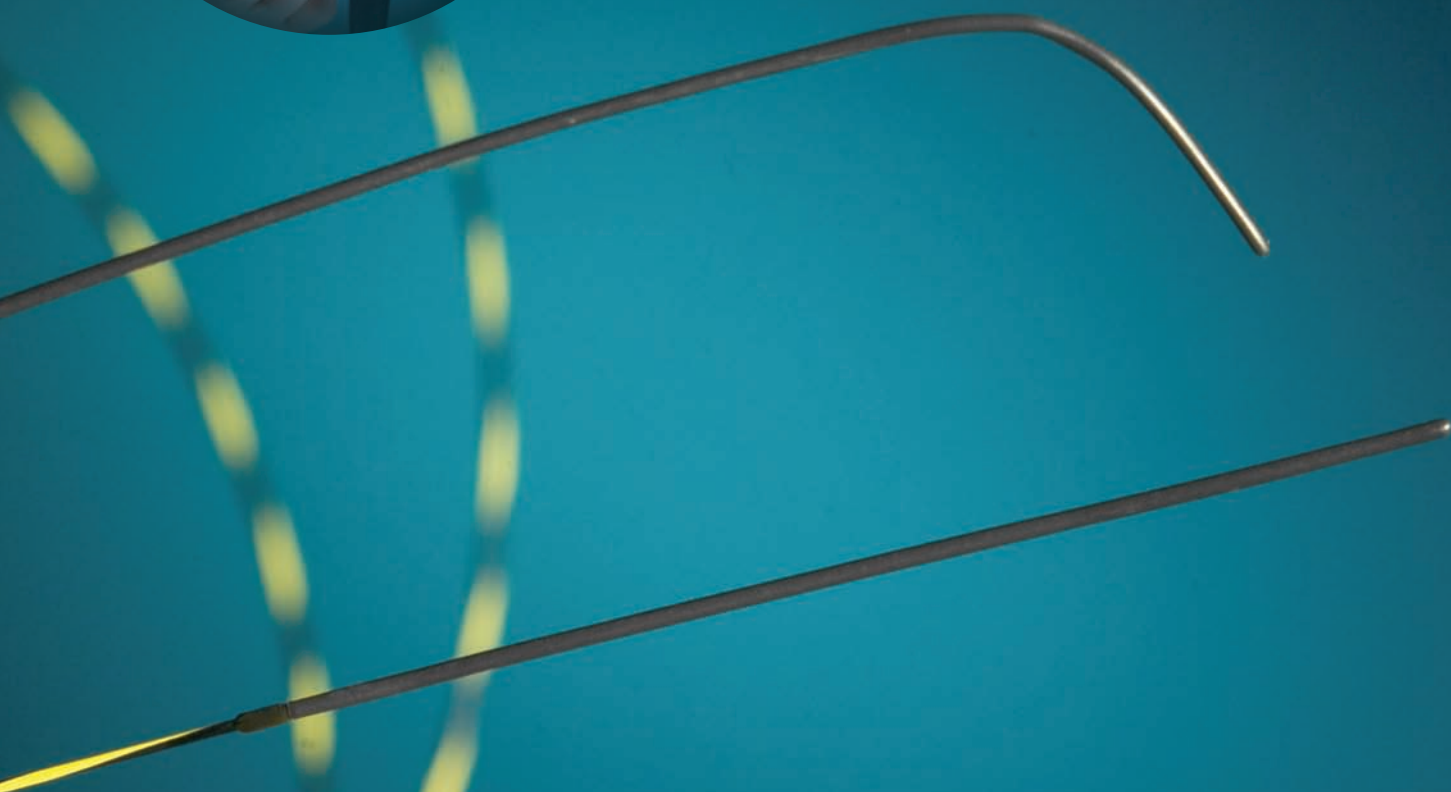

RX Biliary System™

**A difference**  
you can feel.™

# RX Biliary System™

A difference you can feel.™

The endoscope shows you the way, but it's your sense of touch that gets you where you need to go. The devices that make up the RX Biliary System have been developed to provide access and control with a sensation of device responsiveness that truly is... *a difference you can feel.*


Over the years, Boston Scientific has earned the confidence and allegiance of customers around the world who have adopted the RX Biliary System for everyday use. Our legacy of success in developing high performance RX biliary devices has enabled us to change the landscape of ERCP procedures. By innovating and refining products, we've put control and access right at your fingertips.

## TWO DISTINCT TIPS ON A SINGLE GUIDEWIRE

*Guidewire tip technology designed to improve procedural efficiency*

- ◆ 10cm hydrophilic Dream Tip™ End is designed for soft, smooth access and excellent tactile feel. Switch to the 5cm hydrophilic tip to meet procedural demands for firmer tip
- ◆ 5cm tip similar in performance to the Jagwire™ Guidewire
- ◆ 260cm and 450cm guidewire lengths available to meet physician preference for short-wire or traditional over-the-wire technique
- ◆ Guidewire may be left in place during sphincterotomy
- ◆ Choice of angled or straight Dream Tip End

Hyd

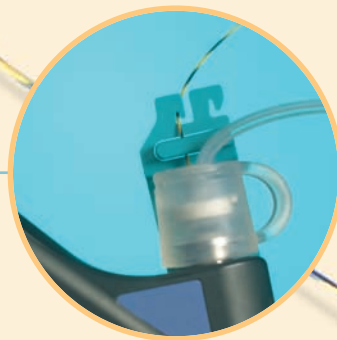

## RX LOCKING DEVICE

*Guidewire locks in place to maintain access*

- ◆ Guidewire may be locked and unlocked at any time during the procedure
- ◆ Facilitates device exchange, contrast injection and biliary therapy



## TACTILE RESPONSIVENESS


*Side-loading Torque Device designed for effective guidewire tip manipulation*

- ◆ Engineered to facilitate deep cannulation of the common bile duct
- ◆ Endoglide™ Coating promotes smooth tracking and enhances tactile sensation
- ◆ Kink-resistant, nitinol shaft engineered for wire control and manipulation capabilities

# Hydra Jagwire™ Guidewire



Choice of Angled  
or Straight Tip

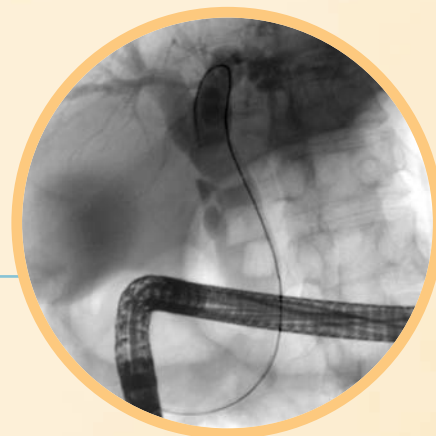

## VISUALISATION

*Hydra Jagwire Guidewire is designed to form an alpha loop to facilitate duct access*

- ◆ Striped movement markers designed to indicate wire positioning during wire advancement and device exchange
- ◆ Radiopaque markers at 10cm and 15cm facilitate accurate stricture measurement
- ◆ Tungsten-filled tips aid fluoroscopic visualisation to help confirm access during device exchange

*Images courtesy of Thomas Kowalski, MD, Thomas Jefferson University Hospital, Philadelphia, PA*

# Hydra Jagwire™

## Guidewire

### HYDRA JAGWIRE GUIDEWIRE

Order Number	Description	Length (cm)	Dream Tip Style	Shaft	Diameter
M00556001	Hydra Jagwire ST	260	Straight		.032"/ 0.89mm
M00556011	Hydra Jagwire ANG	260	Angled		.032"/ 0.89mm
M00556021	Hydra Jagwire ST SS	260	Straight	Stiff	.032"/ 0.89mm
M00556031	Hydra Jagwire ANG SS	260	Angled	Stiff	.032"/ 0.89mm
M00556041	Hydra Jagwire ST	450	Straight		.032"/ 0.89mm
M00556051	Hydra Jagwire ANG	450	Angled		.032"/ 0.89mm
M00556061	Hydra Jagwire ST SS	450	Straight	Stiff	.032"/ 0.89mm
M00556071	Hydra Jagwire ANG SS	450	Angled	Stiff	.032"/ 0.89mm

Packaged 2 per box

#### European Distribution Centre – The Netherlands

T: +31 45 54 67 700 F: +31 45 54 67 800

#### European Headquarters – Paris

T: +33 1 57 66 80 00 F: +33 1 57 66 84 99

#### Argentina (Freephone)

T: 0800 555 2678 F: +5411 4896 8517

#### Australia/New Zealand

T: +61 2 8063 8100 F: +61 2 9330 1404

#### Austria

T: +43 1 60 810 F: +43 1 60 810 60

#### Belgium (Freephone)

T: 0800 94 494 F: 0800 93 343

#### Brazil

T: +55 11 5502 8500 F: +55 11 5103 2212

#### Canada

T: +1 888 359 9691 F: +1 888 575 7396

#### Chile

T: +562 445 4904 F: +562 445 4915

#### China – Beijing

T: +86 10 8525 1588 F: +86 10 8525 1566

#### China – Guangzhou

T: +86 20 8767 9791 F: +86 20 8767 9789

#### China – Shanghai

T: +86 21 6141 5959 F: +86 21 6141 5900

#### Colombia

T: +57 1 629 5045 F: +57 1 612 4761

#### Czech Republic

T: +420 296 331 901 F: +420 296 331 935

#### Denmark (Freephone)

T: 80 30 80 02 F: 80 30 80 05

#### Eire (Freephone)

T: 1800 882 969 F: 1800 882 968

#### Finland

T: +358 20 762 8882 F: +358 20 762 8883

#### France

T: +33 1 39 30 49 00 F: +33 1 39 30 49 01

#### Germany (Freephone)

T: 08000 723300 F: 08000 723319

#### Greece

T: +30 210 9542 300 F: +30 210 9542 310

#### Hong Kong

T: +852 2960 7100 F: +852 2563 5276

#### Hungary

T: +36 1 456 30 40 F: +36 1 456 30 41

#### India – Bangalore

T: +91 80 2212 4928/9 F: +91 80 2207 5153

#### India – Chennai

T: +91 44 2220 1879 F: +91 44 2220 0716

#### India – Delhi

T: +91 11 4243 2222 F: +91 11 2610 0808

#### India – Mumbai

T: +91 22 4030 9165 F: +91 22 4040 9199

#### Italy

T: +39 010 60 60 1 F: +39 010 60 60 200

#### Korea

T: +82 2 3476 2121 F: +82 2 3476 1776

#### Mexico

T: +52 55 5687 63 90 F: +52 55 5687 62 28

#### Middle East/Gulf/North Africa

T: +961 1 805 410 F: +961 1 805 445

#### The Netherlands

T: +31 30 602 55 44 F: +31 30 602 55 05

#### Norway (Freephone)

T: 800 104 04 F: 800 101 90

#### Poland

T: +48 22 435 14 14 F: +48 22 435 14 10

#### Portugal

T: +351 1 381 25 40 F: +351 21 381 25 58

#### South Africa

T: +27 11 840 8600 F: +27 11 463 6077

#### South East Asia – Malaysia

T: +60 3 2283 3813 F: +60 3 2284 3813

#### South East Asia – Philippines

T: +63 2 687 6994 F: +63 2 687 3047

#### South East Asia – Singapore

T: +65 6418 8888 F: +65 6418 8899

#### South East Asia – Thailand

T: +66 2 6380 100 F: +66 2 6380 400

#### South East Europe

T: +30 210 9542 300 F: +30 210 9542 310

#### Spain – Barcelona

T: +34 93 444 72 00 F: +34 93 405 90 45

#### Spain – Madrid

T: +34 90 11 12 15 F: +34 91 319 50 03

#### Sweden

T: +46 42 25 69 00 F: +46 42 25 69 69

#### Switzerland

T: +41 32 626 57 00 F: +41 32 626 57 01

#### Taiwan

T: +886 2 2747 7278 F: +886 2 2747 7270

#### Turkey – Istanbul

T: +90 216 464 36 66 F: +90 216 464 36 67

#### UK (Freephone)

T: 0844 800 4512 F: 0844 800 4513

#### Uruguay

T: +59 82 900 6212 F: +59 82 900 6212

#### Venezuela

T: +58 212 959 6275 F: +58 212 959 5328

**Boston  
Scientific**

*Delivering what's next.™*

[www.bostonscientific.com](http://www.bostonscientific.com)

[www.bostonscientific-international.com](http://www.bostonscientific-international.com)