


Agile™

Esophageal Fully Covered, Removable & Partially Covered Stent System

Featuring a Through-the-Scope
Delivery System for 14mm and 18mm
diameter stents

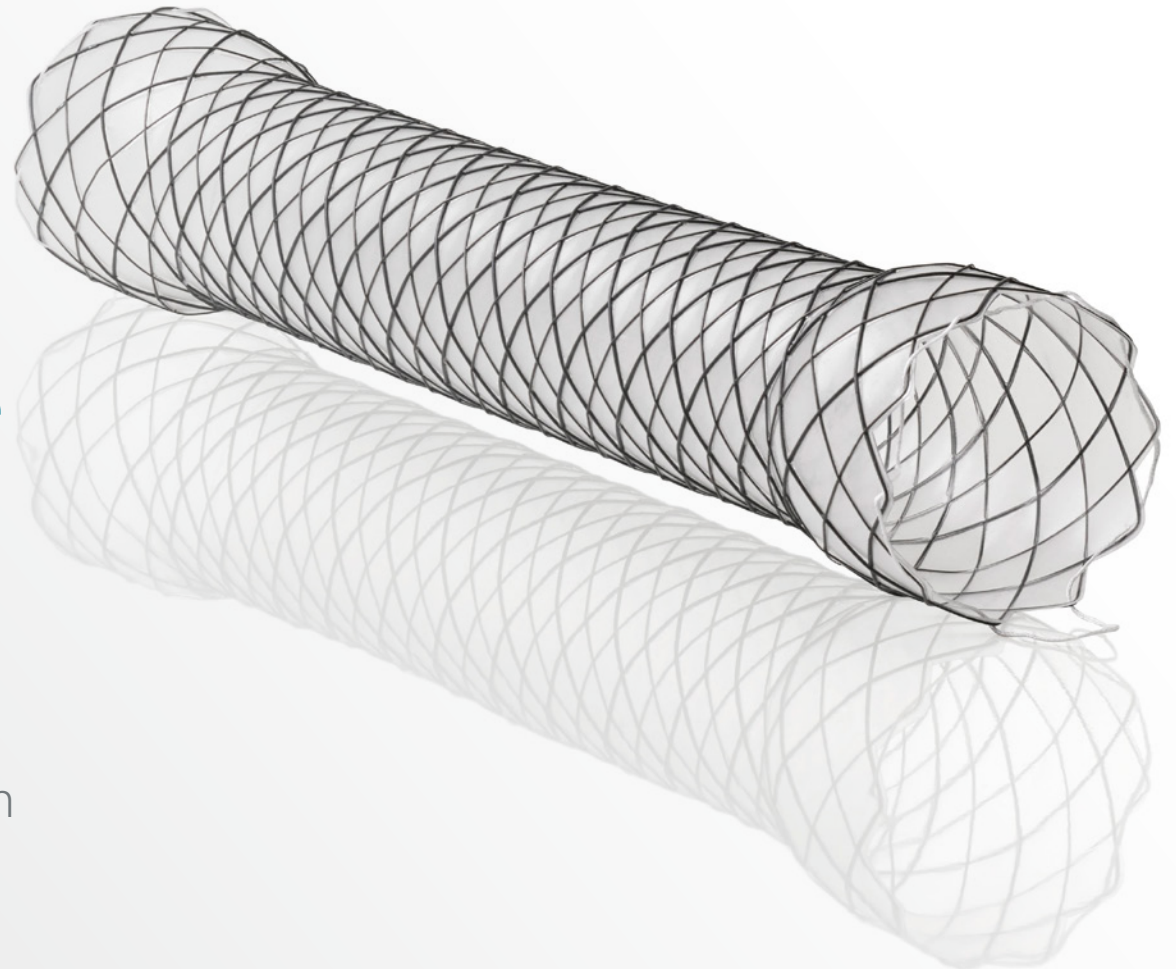


About Agile™

The Agile Esophageal Stent System offers a flexible and conformable self-expanding metal stent on a through the scope delivery system for 14mm and 18mm diameter stents which enables direct visualisation during deployment and may help minimise reliance on fluoroscopy. The 23mm is provided on the trusted WallFlex™ Esophageal over-the-wire delivery system. The Agile portfolio is offered in fully covered and partially covered designs.

- Offering new sizing options including a 60mm length and 14mm diameter designed to meet a variety of clinical needs, all size configurations are offered in partially covered and fully covered removable options.
- Braided dual gauge nitinol wire provides an implant that is more flexible and conformable.*
- Atraumatic looped ends to help minimise risk of trauma.
- Retrieval suture on both ends of stent.
- Tissue In-Growth Prevention – the silicone covering extends the entire length of the Fully Covered Agile Esophageal Stent and is designed to address tumour ingrowth.

*Compared to WallFlex™ Esophageal 18mm and 23mm stent diameters. WallFlex Esophageal is not available in a 14mm stent diameter for direct comparison.

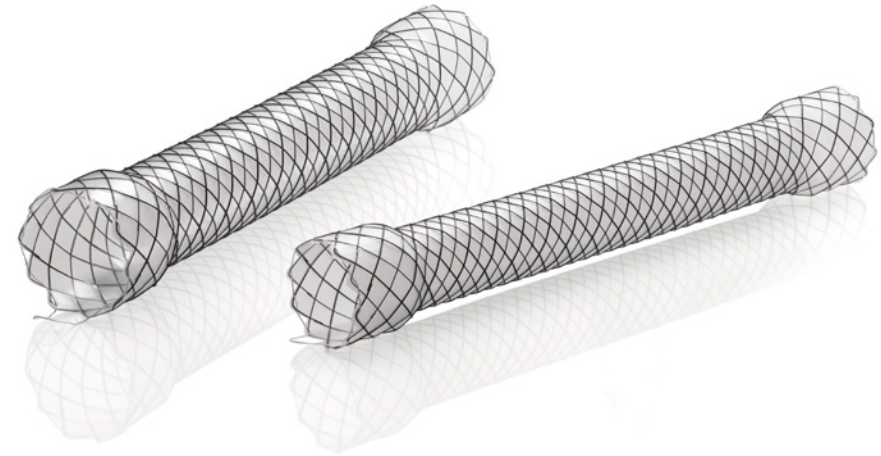

Delivery System

Through-the-Scope design

- For 14mm and 18mm stents, the through-the-scope delivery system provides direct visualization to aid physicians in stent placement accuracy and help minimize fluoroscopic dependence.
- Designed to offer flexibility and trackability to enhance physician control.

Low profile delivery system


Compatible with a 3.7mm working channel, designed to navigate through tight strictures.

5 radiopaque markers

Markers signal the ends of the stent, the midpoint of the stent, and the point at which the stent can no longer be recaptured, aiding in deployment. The fifth marker at the leading end of the exterior tube indicates how far the stent has been deployed.

Reconstrainability

In the event that the stent is not exactly where a physician intends it to be, the stent can be reconstrained up to 2x and redeployed.

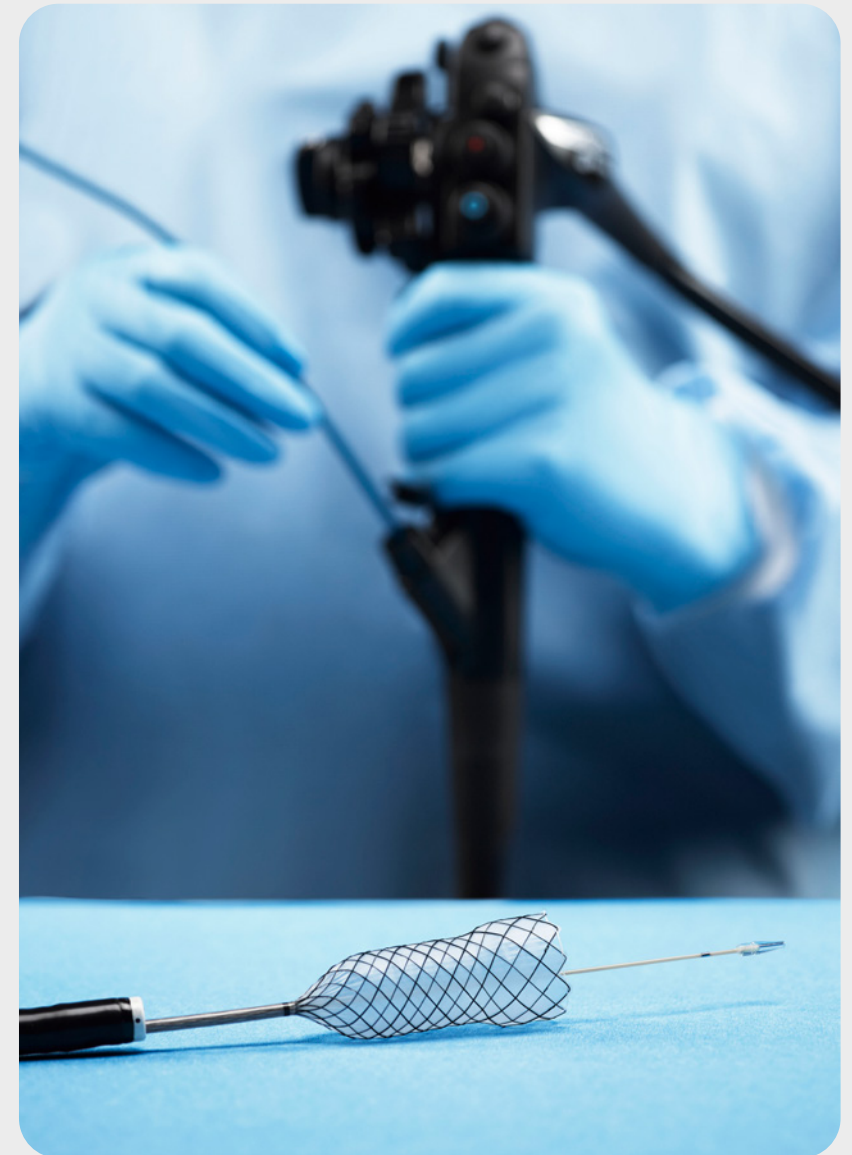

Ordering Information

The Agile™ Esophageal Partially Covered Stent System is intended for maintaining esophageal luminal patency in esophageal strictures caused by intrinsic and/or extrinsic malignant tumours and occlusion of concurrent esophageal fistulas.

Agile Esophageal Stent System – Partially Covered					
Order Number	GTIN	Stent O.D. (mm)	Proximal/Distal Flares O.D. (mm)	Stent Length (cm)	Delivery System Diameter (Fr)
M00517200	08714729972907	14	19	6.2	10.5
M00517210	08714729972914	14	19	10.2	10.5
M00517220	08714729972921	14	19	11.9	10.5
M00517230	08714729972938	14	19	14.8	10.5
M00517240	08714729972945	18	23	5.9	10.5
M00517250	08714729972952	18	23	9.7	10.5
M00517260	08714729972969	18	23	11.9	10.5
M00517270	08714729972976	18	23	14.9	10.5

Agile Esophageal Stent System – OTW Partially Covered					
M00517710	00191506006587	23	28	6.2	18.5
M00517720	00191506006594	23	28	10.1	18.5
M00517730	00191506006600	23	28	12.0	18.5
M00517740	00191506006617	23	28	15.0	18.5



Ordering Information

The Agile™ Esophageal Fully Covered Removable (RMV) Stent System is intended for maintaining esophageal luminal patency in esophageal strictures caused by intrinsic and/or extrinsic malignant tumours, occlusion of concurrent esophageal fistulas and treating refractory benign esophageal strictures.

Agile Esophageal Stent System – Fully Covered – Removable					
Order Number	GTIN	Stent O.D. (mm)	Proximal/Distal Flares O.D. (mm)	Stent Length (cm)	Delivery System Diameter (Fr)
M00517580	08714729973133	14	19	6.2	10.5
M00517590	08714729973140	14	19	10.2	10.5
M00517600	08714729973157	14	19	11.9	10.5
M00517610	08714729973164	14	19	14.8	10.5
M00517620	08714729973171	18	23	5.9	10.5
M00517630	08714729973188	18	23	9.7	10.5
M00517640	08714729973195	18	23	11.9	10.5
M00517650	08714729973201	18	23	14.9	10.5

Agile Esophageal Stent System – OTW Fully Covered Removable					
M00517790	00191506006662	23	28	6.2	18.5
M00517800	00191506006679	23	28	10.1	18.5
M00517810	00191506006686	23	28	12.0	18.5
M00517820	00191506006693	23	28	15.0	18.5



All photographs owned by Boston Scientific.

All cited trademarks are the property of their respective owners. CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labelling supplied with each device. Products shown for INFORMATION purposes only and may not be approved or for sale in certain countries.

This material not intended for use in France.

Product available in the European Economic Area (EEA) only. Please check availability with your local sales representative or customer service.

**Boston
Scientific**
Advancing science for life™

www.bostonscientific.eu

© 2022 Boston Scientific Corporation or its affiliates. All rights reserved.

ENDO-1093104-AC

CE 0344