

WOVEN™ / WOVEN FLEXIE™

Fixed Curve Diagnostic Catheters

Design for reliability and predictability

Technical Information

Description	Specifications
Shaft Diameter	4F (1.33 mm), 5F (1.66 mm), 6F (2 mm) and 7F (2.33 mm)
Electrode Spacing	1 mm, 2 mm, 5 mm, 10 mm, 2,5,2 mm, 5,5,5,282 mm, 10,10,267 mm, 5,5,5,152 mm, 1,10,1 mm, 2,7,2 mm, 2,9,2 mm, 5,5,380 mm
Usable Length	80 cm, 100 cm, 110 cm, 120 cm, 125 cm, 150 cm
Electrode Material	Pure Pt, Ag connection wire
Tip Options	See charts
Conguration	2e, 3e, 4e, 5e, 6e, 10e, 2e, 24e

Key Features

Platinum electrodes connected to pure silver core conducting wires, used only by the WOVEN™ range to produce the highest resolution electrograms

The WOVEN™ catheter's bioexibility allows the shaft to soften rapidly at body temperature to a predetermined biosoftness

Biosoftness characteristics allows the catheter tip to conform to the endocardial wall for excellent tissue contact


Softening response reduces endocardial bounce, providing greater stability and better signal quality

The primary difference between WOVENFLEXIE™ and the standard WOVEN™ catheter is the variable durometer coating system, which increases tip exibility by 17%

Interweaving 16 strands of pure silver core conducting wire provides excellent conduction with low impedance


WOVEN™

Fixed Curve Diagnostic Catheter

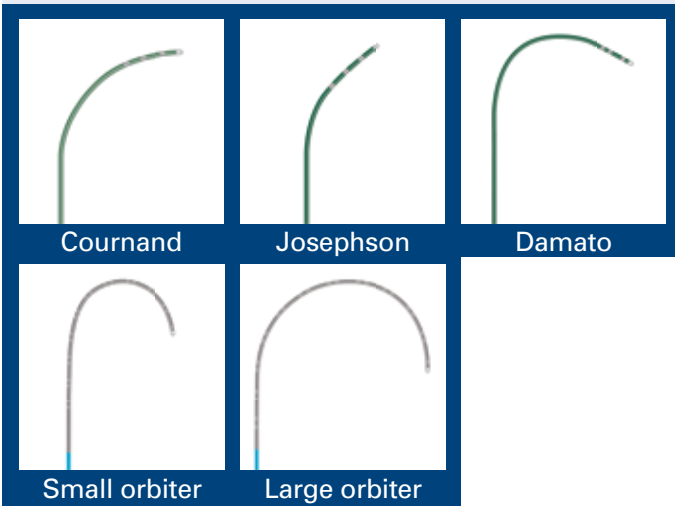

WOVEN FLEXIE™

Fixed Curve Diagnostic Catheter

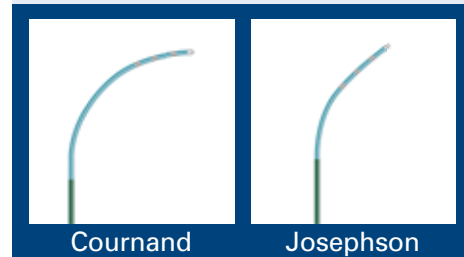

Curve Options

Curve Options

WOVEN™ / WOVEN FLEXIE™

Fixed Curve Diagnostic Catheters

Ordering Information

WOVEN™ BIPOLAR

Catalog number	French size	N° of electrodes	Curve style	Electrode spacing	Usable length	Cables
M0041910010	6F (2 mm)	2	Cournand	10 mm	125 cm	M004200088P0

WOVEN™ TRIPOLAR

Catalog number	French size	N° of electrodes	Curve style	Electrode spacing	Usable length	Cables
M004002853S0	7F (2.33 mm)	3	Damato	10 mm	110 cm	See Chart Below*

WOVEN™ QUADPOLAR

Catalog number	French size	N° of electrodes	Curve style	Electrode spacing	Usable length	Cables
M0042001500	4F (1.33 mm)	4	Cournand	2,5,2 mm	80 cm	M004008568P0
M0042001510	4F (1.33 mm)	4	Josephson	2,5,2 mm	80 cm	M004008568P0
M0044FMC005440	4F (1.33 mm)	4	Josephson	5 mm	110 cm	M004200088P0
M0042001530	4F (1.33 mm)	4	Josephson	1 mm	80 cm	M004008568P0
M0044FMC003030	4F (1.33 mm)	4	Josephson	2,5,2 mm	100 cm	M004008568P0
M004200056E0	5F (1.66 mm)	4	Cournand	10 mm	110 cm	M004200088P0
M004200580E0	5F (1.66 mm)	4	Cournand	5 mm	110 cm	M00420088P0
M0042005810	5F (1.66 mm)	4	Cournand	2,5,2 mm	110 cm	M004008568P0
M0042000600	5F (1.66 mm)	4	Josephson	5 mm	110 cm	M004008568P0
M004200060E0	5F (1.66 mm)	4	Josephson	5 mm	110 cm	M004200088P0
M004200064P0	5F (1.66 mm)	4	Josephson	2,5,2 mm	110 cm	Tails connector only. No cable required.
M0042014000	5F (1.66 mm)	4	Damato	10 mm	110 cm	M004008568P0
M0042005820	6F (2 mm)	4	Cournand	2,5,2 mm	110 cm	M004008568P0
M0042005830	6F (2 mm)	4	Cournand	2,5,2 mm	125 cm	M004200088P0
M004006798P0	6F (2 mm)	4	Josephson	5 mm	120 cm	M00420088P0
M0042005840	6F (2 mm)	4	Josephson	2,5,2 mm	125 cm	M004200088P0
M0042005910	6F (2 mm)	4	Josephson	2,5,2 mm	125 cm	M004200088P0
M0046FMC004680	6F (2 mm)	4	Josephson	5 mm	120 cm	M004200088P0

WOVEN™ QUADPOLAR WITH SURELINK CONNECTORS

Catalog number	French size	N° of electrodes	Curve style	Electrode spacing	Usable length	Cables
M004006245S0	6F (2 mm)	4	Cournand	10 mm	110 cm	See Chart Below*
M004200066S0	6F (2 mm)	4	Cournand	5 mm	125 cm	See Chart Below*
M004200201S0	6F (2 mm)	4	Josephson	10 mm	125 cm	See Chart Below*
M004200624S0	6F (2 mm)	4	Josephson	5 mm	120 cm	See Chart Below*
M0046FMC006770	6F (2 mm)	4	Josephson	5 mm	120 cm	See Chart Below*
M0046FSC0007S0	6F (2 mm)	4	Josephson	5 mm	120 cm	See Chart Below*
M004200719S0	6F (2 mm)	4	Damato	10 mm	110 cm	See Chart Below*

WOVEN™ / WOVEN FLEXIE™

Fixed Curve Diagnostic Catheters

WOVEN™ PENTAPOLAR

Catalog number	French size	N° of electrodes	Curve style	Electrode spacing	Usable length	Cables
M0042014010	6F (2 mm)	5	Josephson	5,5,5,282 mm	120 cm	M004200089P0
M0042005940	6F (2 mm)	5	Josephson	5,5,5,152 mm	120 cm	M004200089P0

WOVEN™ HEXAPOLAR

Catalog number	French size	N° of electrodes	Curve style	Electrode spacing	Usable length	Cables
M0042002040	5F (1.66 mm)	6	Cournand	5 mm	110 cm	M004006590P0
M0046FMC008210	6F (2 mm)	6	Cournand	5 mm	125 cm	M004560004A0

WOVEN™ DECAPOLAR

Catalog number	French size	N° of electrodes	Curve style	Electrode spacing	Usable length	Cables
M0042003550	6F (2 mm)	10	Cournand	5 mm	125 cm	M004006590P0
M0042002110	6F (2 mm)	10	Cournand	2,5,2 mm	125 cm	M004006590P0
M0042002120	6F (2 mm)	10	Cournand	2 mm	125 cm	M004006590P0
M0046FMC007350	6F (2 mm)	10	Cournand	5 mm	125 cm	M004560004A0
M0046FMC007020	6F (2 mm)	10	Cournand	2,5,2 mm	125 cm	M004560004A0
M0046FMC002420	6F (2 mm)	10	Josephson	1,10,1 mm	125 cm	M004006590P0
M0046FMC008090	6F (2 mm)	10	Josephson	1,10,1 mm	125 cm	M004560004A0

WOVEN™ 12-ELECTRODES

Catalog number	French size	N° of electrodes	Curve style	Electrode spacing	Usable length	Cables
M0042006930	6F (2 mm)	12	Cournand	10 mm	125 cm	M004200089P0 x 2 cables
M0046FMT0204P0	6F (2 mm)	12	Cournand	10 mm	125 cm	Tails connector only. No cable required.

WOVEN™ / WOVEN FLEXIE™

Fixed Curve Diagnostic Catheters

WOVEN™ 24-ELECTRODES

Catalog number	French size	N° of electrodes	Curve style	Electrode spacing	Usable length	Cables
M0042424030	6F (2 mm)	24	Orbiter Small	2,7,2 mm	110 cm	M004200774P0 x 3 cables
M0042424050	6F (2 mm)	24	Orbiter Small	2,9,2 mm	110 cm	M004200774P0 x 3 cables
M0042424040	6F (2 mm)	24	Orbiter Large	2,7,2 mm	110 cm	M004200774P0 x 3 cables
M0042424060	6F (2 mm)	24	Orbiter Large	2,9,2 mm	110 cm	M004200774P0 x 3 cables
M0046FMC007980	6F (2 mm)	24	Orbiter Large	2,9,2 mm	100 cm	M004560004A0 x 3 cables

WOVEN™ FLEXIE QUADPOLAR AND PENTAPOLAR

Catalog number	French size	N° of electrodes	Curve style	Electrode spacing	Usable length	Cables
M0042005960	5F (1.66 mm)	4	Cournand	5 mm	120 cm	M004200088P0
M0042005970	5F (1.66 mm)	4	Josephson	5 mm	120 cm	M004200088P0
M0042005980	5F (1.66 mm)	4	Josephson	2,5,2 mm	120 cm	M004200088P0
M0045FMC007920	5F (1.66 mm)	5	Josephson	5,5,380 mm	110 cm	M004200088P0
M0045002020	6F (2 mm)	4	Cournand	5 mm	120 cm	M004200088P0
M0045002030	6F (2 mm)	4	Josephson	5 mm	120 cm	M004200088P0

Surelink Cable	Description	Cable length
M004560002A0	Protected 4-pin cable	120 cm
M004560002BA0	Protected 4-pin (Blue)	120 cm
M004560002RA0	Protected 4-pin (Red)	120 cm
M004560002YA0	Protected 4-pin (Yellow)	120 cm

**Boston
Scientific**

Advancing science for life™

www.bostonscientific.eu

All cited trademarks are the property of their respective owners. CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labeling supplied with each device. Information for the use only in countries with applicable health authority product registrations. Information not intended for distribution in France.

EP-401203-AA JUN2016 Printed in Germany by medicalvision.

© 2017 Boston Scientific Corporation or its affiliates. All rights reserved.
DINEPxxxxEA