

WOVEN/WOVENFLEXIE™


FIXED CURVE DIAGNOSTIC CATHETERS

One of The first diagnostic EP catheters and a trusted choice by Electrophysiologists for over 25 years.

Designed for reliability and predictability.


Technical Information

Description	Specifications
Shaft Diameter	4F, 5F, 6F and 7F
Usable Length	65cm, 80cm, 90cm, 100cm, 110cm, 115cm, 120cm, 125cm
Electrode Material	Platinum
Curve Size	See Charts
Tip Options	Woven-Fixed Tip WovenFlexie-Soft Tip
Configuration	Quadripolar/Hexapolar/ Pentapolar/Decapolar

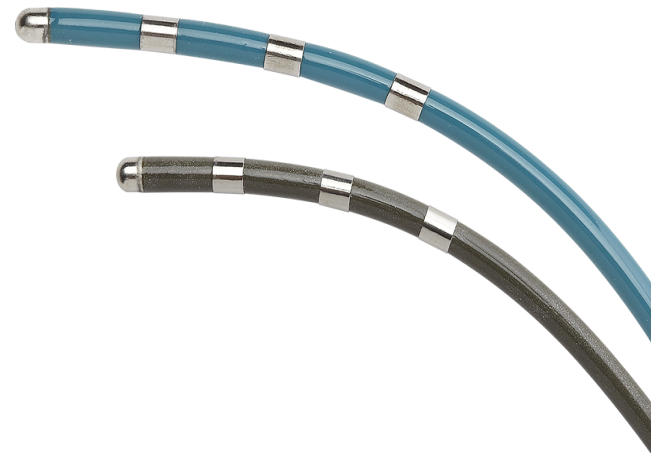


Curve Options

Woven Fixed Curve Diagnostic Catheter

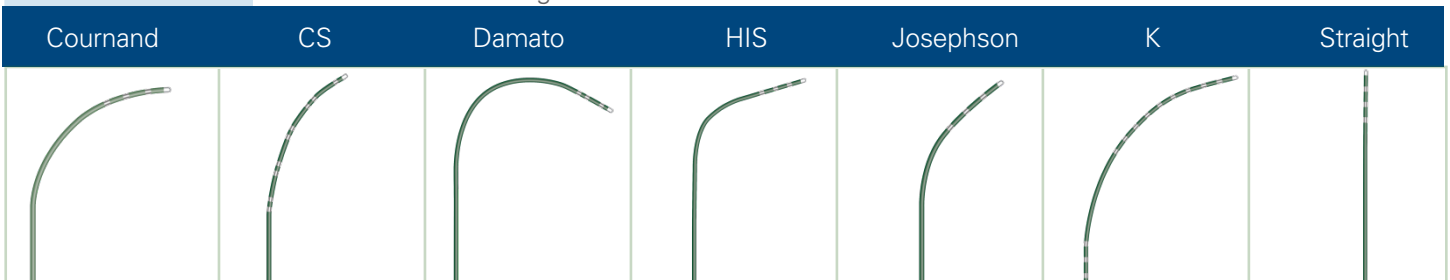

Curve Options

WovenFlexie Fixed Curve Diagnostic Catheter

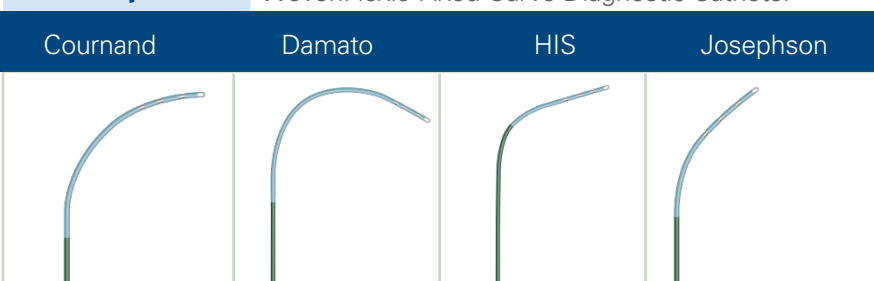

Referenced catheter configurations are illustrative representations only and may not reflect actual performance.

WOVEN/WOVENFLEXIE™

FIXED CURVE DIAGNOSTIC CATHETERS

Ordering Information

Woven Fixed Curve Diagnostic Catheter - 4F

Catheter Model No.	Curve Size	Electrode Configuration	Electrode Spacing	Usable Length	Cable Model
M004 200151 0	Josephson	Quadripolar	2,5,2mm	80cm	M004 008568P
M004 4FMC00544 0	Josephson	Quadripolar	5mm	100cm	M004 200088P

Ordering Information

Woven Fixed Curve Diagnostic Catheter - 5F

Catheter Model No.	Curve Size	Electrode Configuration	Electrode Spacing	Usable Length	Cable Model
M004 5FMC00496 0	Josephson	Quadripolar	2,5,2mm	100cm	M004 200088P
M004 200060E 0	Josephson	Quadripolar	5mm	110cm	M004 200088P
M004 200060 0	Josephson	Quadripolar	5mm	110cm	M004 008568P
M004 200577E 0	Josephson	Quadripolar	5mm	120cm	M004 200088P
M004 200580E 0	Cournand	Quadripolar	5mm	110cm	M004 200088P
M004 200056E 0	Cournand	Quadripolar	10mm	110cm	M004 200088P
M004 200581 0	Cournand	Quadripolar	2,5,2mm	110cm	M004 008568P

Ordering Information

Woven Fixed Curve Diagnostic Catheter - 6F

Catheter Model No.	Curve Size	Electrode Configuration	Electrode Spacing	Usable Length	Cable Model
M004 191001 0	Cournand	Bipolar	10mm	125cm	M004 200088P
M004 200584 0	Josephson	Quadripolar	2,5,2mm	125cm	M004 200088P
M004 200624S 0	Josephson	Quadripolar	5mm	120cm	M004 560002A
M004 6FSC00007 0	Josephson	Quadripolar	5mm	120cm	M004 200088P
M004 6FSC00075 0	Josephson	Quadripolar	5mm	120cm	M004 560002A
M004 6FMC00468 0	Josephson	Quadripolar	5mm	120cm	M004 200088P
M004 6FMC006770	Josephson	Quadripolar	5mm	120cm	M004 560002A
M004 200066S 0	Cournand	Quadripolar	5mm	125cm	M004 560002A
M004 200201S 0	Josephson	Quadripolar	10mm	125cm	M004 560002A
M004 200583 0	Cournand	Quadripolar	2,5,2mm	125cm	M004 200088P
M004 200574 0	Cournand	Quadripolar	5mm	125cm	M004 200088P
M004 006245S 0	Cournand	Quadripolar	10mm	110cm	M004 560002A
M004 2001123 0	Cournand	Quadripolar	10mm	110cm	M004 200088P
M004 200586 0	Damato	Quadripolar	2,5,2mm	125cm	M004 200088P
M004 200594 0	Josephson	Pentapolar	5,5,5,152mm	120cm	M004 200089P
M004 201401 0	Josephson	Pentapolar	5,5,5,282mm	120cm	M004 200089P
M004 6FMC00821 0	Cournand	Hexapolar	5mm	125cm	M004 560004A
M004 6FMC00823 0	Damato	Hexapolar	2,10,2mm	90cm	M004 560004A

WOVEN/WOVENFLEXIE™

FIXED CURVE DIAGNOSTIC CATHETERS

Ordering Information

Woven Fixed Curve Diagnostic Catheter - 6F

Catheter Model No.	Curve Size	Electrode Configuration	Electrode Spacing	Usable Length	Cable Model
M004 200211 0	Cournand	Decapolar	2,5,2mm	125cm	M004 006590P
M004 200212 0	Cournand	Decapolar	2mm	125cm	M004 006590P
M004 200355 0	Cournand	Decapolar	5mm	125cm	M004 006590P
M004 200578 0	Damato	Decapolar	2,5,2mm	110cm	M004 006590P
M004 6FMC00735 0	Cournand	Decapolar	5mm	125cm	M004 560004A
M004 6FMC00809 0	Josephson	Decapolar	1,10,1mm	125cm	M004 560004A
M004 200693 0	Cournand	12E	10mm	125cm	Requires 2x200774P
M004 6FMC00798 0	Orbiter Large	24E	2,9,2mm	100cm	2ea. 560001A 1ea. 560002A

Ordering Information

Woven Fixed Curve Diagnostic Catheter - 7F

Catheter Model No.	Curve Size	Electrode Configuration	Electrode Spacing	Usable Length	Cable Model
M004 7FMC00675 0	Josephson	Decapolar	5mm	115cm	M004 560004A

Ordering Information

WovenFlexie Fixed Curve Diagnostic Catheter -5F

Catheter Model No.	Curve Size	Electrode Configuration	Electrode Spacing	Usable Length	Cable Model
M004 200597 0	Josephson	Quadripolar	5mm	120cm	M004 200088P
M004 5FMC00711 0	Josephson	Quadripolar	5mm	125cm	M004 560002A
M004 5FMC00792 0	Josephson	Hexapolar	5,5,380mm	110cm	M004 200089P

Ordering Information

WovenFlexie Fixed Curve Diagnostic Catheter - 6F

Catheter Model No.	Curve Size	Electrode Configuration	Electrode Spacing	Usable Length	Cable Model
M004 500203 0	Josephson	Quadripolar	5mm	120cm	M004 200088P
M004 500202 0	Cournand	Quadripolar	5mm	120cm	M004 200088P

WOVEN AND WOVENFLEXIE™ FIXED CURVE DIAGNOSTIC CATHETERS

INTENDED USE: BARD® Electrophysiology's fixed curve diagnostic electrode catheters are intended for temporary intracardiac sensing, recording, stimulation and temporary pacing during the evaluation of cardiac arrhythmias.

CONTRAINDICATIONS: The transseptal approach is contraindicated in patients with left atrial thrombus or myxoma, or interatrial baffle patch. The retrograde transaortic approach is contraindicated in patients with aortic valve replacement.

WARNINGS: The risks of using electrophysiology catheters include those risks related to heart catheterization, such as thromboembolism, perforation, tamponade, and infection. The induction of an unintended arrhythmia is a known complication. Do not use excessive force to advance or withdraw the catheter when resistance is encountered. If using an open lumen catheter remove any guidewire/stylette prior to electrical stimulation.

PRECAUTIONS: Excessive bending, torquing, or kinking of the electrode catheter may cause damage to the catheter, including damage to the internal wires. Use only sterile saline or water to wipe this catheter. 91099783 (REV AA)

CAUTION: Federal law (USA) restricts this device to sale by or on the order of a physician. Rx only. Prior to use, please see the complete "Directions for Use" for more information on Indications, Contraindications, Warnings, Precautions, Adverse Events, and Operator's Instructions.



Rhythm Management
300 Boston Scientific Way
Marlborough, MA 01752-1234
www.bostonscientific.com

Medical Professionals:
1.800.CARDIAC (227.3422)
Customer Service:
1.888.272.1001

© 2018 Boston Scientific Corporation
or its affiliates. All rights reserved.

EP-349605-AB