


WallFlex™ Colonic and Duodenal Soft Stent System with Anchor Lock Delivery System

**Boston
Scientific**
Advancing science for life™

Built on the trusted WallFlex™ platform, the WallFlex Colonic and Duodenal Soft Stent Systems are designed to be a more flexible option to relieve malignant obstructions of the colon and duodenum.

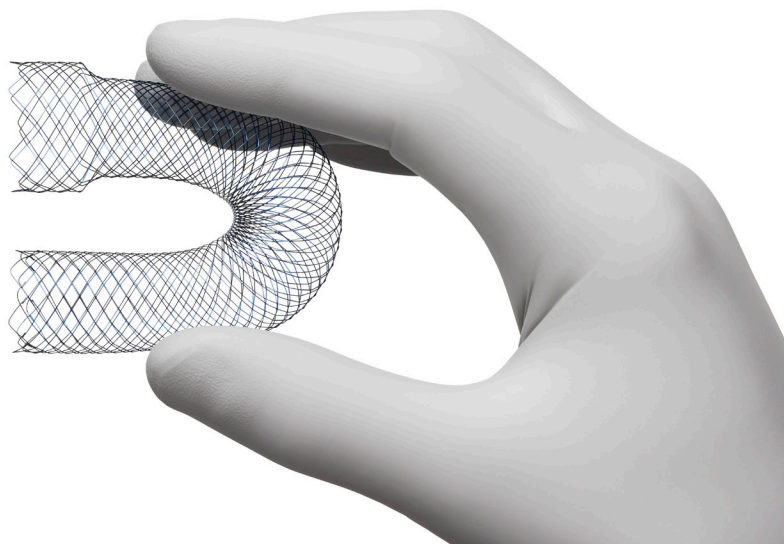


Softer Stent

Designed with a lower axial force which may better conform to torturous anatomy such as the sigmoid, flexures, and duodenal bulb.

Delivery System

New low profile 9Fr delivery system, allowing for placement via a 3.2mm channel scope, which may assist in placement in tighter strictures.

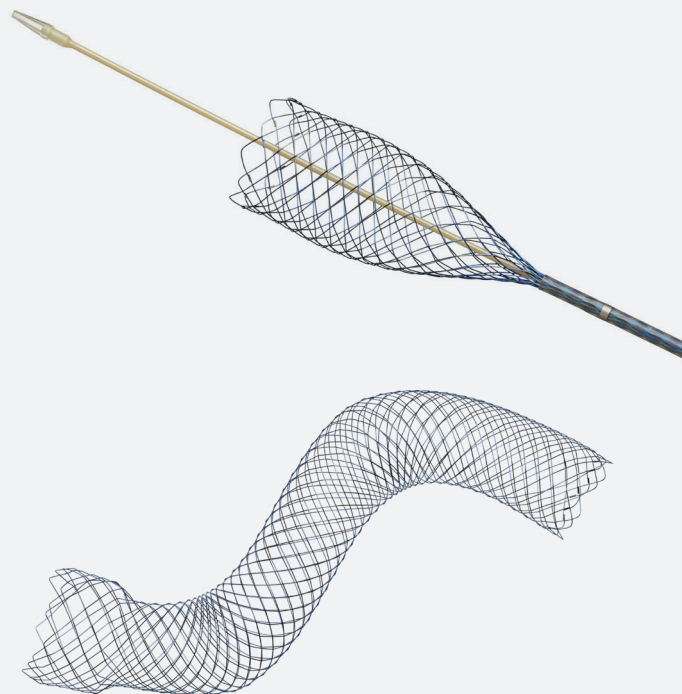

Flexibility

With lower radial and axial forces*, WallFlex Colonic and Duodenal Soft Stents are designed for increased flexibility.

Comprehensive Metal Stent Portfolio

Boston Scientific's robust line of esophageal, biliary and enteral metal stents provides physicians the tools needed to help them effectively manage strictures. WallFlex Enteral stents have been studied in more than 100 clinical publications since 2004.

*as compared to the original WallFlex Colonic and Duodenal Stent System

WallFlex™ Colonic and Duodenal Soft Stent System with Anchor Lock Delivery System

The WallFlex Colonic Soft Stent is indicated for the palliative treatment of colonic strictures caused by malignant neoplasm and to relieve large bowel obstruction prior to colectomy in patients with malignant strictures. The WallFlex Enteral Duodenal Stent Soft is indicated for palliative treatment of gastroduodenal obstructions produced by malignant neoplasms.

Ordering Information

WallFlex Colonic Soft Stent System with Anchor Lock Delivery System

Order Number	Diameter (mm) Flare/ Body	Stent Length (cm)	Working Length (cm)	Delivery System Diameter/Working Channel (F/mm)	Guidewire Diameter (inches/mm)
M00555760	27 / 22	6	230	9 / 3.2	0.035 / 0.89
M00555770	27 / 22	9	230	9 / 3.2	0.035 / 0.89
M00555780	27 / 22	12	230	9 / 3.2	0.035 / 0.89

WallFlex Duodenal Soft Stent System with Anchor Lock Delivery System

Order Number	Diameter (mm) Flare/ Body	Stent Length (cm)	Working Length (cm)	Delivery System Diameter/Working Channel (F/mm)	Guidewire Diameter (inches/mm)
M00555860	27 / 22	6	230	9 / 3.2	0.035 / 0.89
M00555870	27 / 22	9	230	9 / 3.2	0.035 / 0.89
M00555880	27 / 22	12	230	9 / 3.2	0.035 / 0.89

**Boston
Scientific**
Advancing science for life™

Boston Scientific Corporation
300 Boston Scientific Way
Marlborough, MA 01752

www.bostonscientific.com/gastro

© 2025 Boston Scientific Corporation
or its affiliates. All rights reserved.

ENDO-815206-AB
CE 0344

CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings, and instructions for use can be found in the product labelling supplied with each device or at www.IFU-BSCI.com. Products shown for INFORMATION purposes only and may not be approved or for sale in certain countries. This material not intended for use in France. Rx Only

All photographs taken by Boston Scientific.
All trademarks are the property of their respective owners.