

Interlock™-35 Fibered IDC™ Occlusion System



Standard Length 2D Configurations				
UPN	Description	Diameter (mm)	Length (cm)	Shape
M001363500	Interlock - 35 Coil	3	4	2D
M001363520	Interlock - 35 Coil	4	10	2D
M001363540	Interlock - 35 Coil	6	10	2D
M001363550	Interlock - 35 Coil	6	20	2D
M001363570	Interlock - 35 Coil	8	10	2D
M001363580	Interlock - 35 Coil	8	20	2D
M001363590	Interlock - 35 Coil	8	40	2D
M001363600	Interlock - 35 Coil	10	20	2D
M001363610	Interlock - 35 Coil	10	40	2D
M001363620	Interlock - 35 Coil	12	20	2D
M001363630	Interlock - 35 Coil	12	40	2D
M001363640	Interlock - 35 Coil	15	20	2D
M001363650	Interlock - 35 Coil	15	40	2D
M001363660	Interlock - 35 Coil	18	20	2D
M001363670	Interlock - 35 Coil	18	40	2D



Cube Configurations				
UPN	Description	Diameter (mm)	Length (cm)	Shape
M001363700	Interlock - 35 Coil	4	6	Cube
M001363720	Interlock - 35 Coil	6	10	Cube
M001363730	Interlock - 35 Coil	6	20	Cube
M001363760	Interlock - 35 Coil	8	20	Cube
M001363790	Interlock - 35 Coil	10	25	Cube
M001363800	Interlock - 35 Coil	10	40	Cube
M001363810	Interlock - 35 Coil	15	25	Cube
M001363820	Interlock - 35 Coil	15	40	Cube
M001363830	Interlock - 35 Coil	20	40	Cube



Diamond Configurations				
UPN	Description	Diameter (mm)	Length (cm)	Shape
M001363910	Interlock - 35 Coil	4	4.5	Diamond
M001363920	Interlock - 35 Coil	6	9	Diamond
M001363930	Interlock - 35 Coil	8	14	Diamond

# Interlock™-35 Coil

Fibered IDC™ Occlusion System



The evolution begins –  
**one advancement  
after another**

# Accuracy, Precision and Power

The Interlock™-35 Coil introduced the unique combination of powerful Dacron™ Fibers and precise detachability. Now we redefine deliverability with advancements that are designed to improve the consistency of coil deployment. Technology advances again.

## Improved Deliverability

Side holes on coil introducer allow for hydration during advancement and retraction

Nitinol pusher wire engineered for improved kink resistance

Flushing luer intended to facilitate pre-deployment hydration

## Technology advances again.

### Precision


Interlocking arms for precise placement

- The Interlock-35 Coil features a simple interlocking connection between the pusher wire and the coil.
- While inside the lumen of the catheter, the interlocking connection remains attached.
- When pushed outside of the catheter, the connection detaches.

### Power

Dacron Fibers for rapid occlusion

#### Number of Coils Used Per Procedure

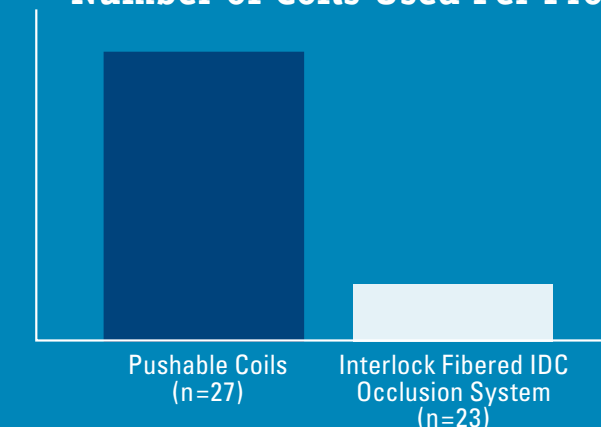

Dudek, et al, Embolization of the Gastroduodenal Artery Before Selective Internal Radiotherapy: A Prospectively Randomized Trial Comparing Standard Pushable Coils with Fibered Interlock Detachable Coils. CVIR, April 14, 2010.

- 83.6% of cases (19 of 23) achieved complete occlusion with one Interlock Fibered IDC Occlusion Coil.
- A shorter procedure time was observed for procedures utilizing detachable coils when compared to pushable coils.
- No coil migration was observed in the detachable coil group, whereas one incidence of migration was observed in the pushable coil group.