

# Swiss LithoClast® Trilogy Lithotripter

**Dual-energy, single-probe system.** The LithoClast Trilogy Lithotripter delivers controlled ultrasonic and ballistic energy concurrently or independently with simultaneous suction through a single probe.

**1. Centralized controls  
with touchscreen interface.**

Designed for an intuitive user experience to customize energy and suction settings.

**2. StoneCatcher.**

A proprietary suction tubing and stone collection system designed for seamless collection of stone fragments for analysis.

**3. Pistol-grip handpiece.**

An ergonomical pistol-grip hand piece designed for prone, supine, and split-leg positioning.



**Smart Technology.**

System automatically regulates maximum energy and suction settings for each probe. **Diagnose and solve common errors through built-in diagnostic and troubleshooting tools.** Track ultrasound performance with an ultrasound efficiency meter during the procedure. **RFID tags in the probe and handpiece store the component's product and usage information.**



## Plug and play system.

Easy set up with minimal training required.

### 4. Integrated pump with suction control.

A peristaltic pump with adjustable settings controls the suction flow rate during energy activation. Suction activation is controlled by the physician.

### 5. Foot switch.


Control energy and suction with a single pedal. A half step activates suction, and a full step triggers suction and energy.

### 6. Suction bags.

Disposable suction collection bags that turn liquid into solid biomedical waste for safe disposal.

### 7. Disposable Probes.

A wide array of single-use probe sizes compatible with most rigid scopes.\* Color coded for quick identification. Offered with a torque wrench for optimum probe tightening

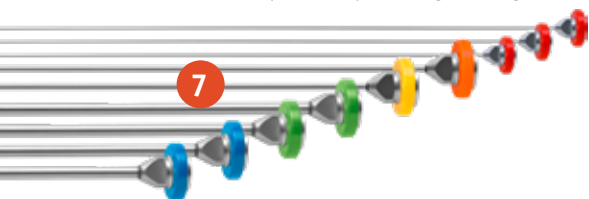

### 8. Mobile cart.\*

Mount the system on a compact cart for portability.

### 9. HDMI/USB ports.\*

Project console screen on secondary monitor for physician convenience. Usage files can be downloaded by a trained technician.



**Boston  
Scientific**

Advancing science for life™

Manufactured by:  
EMS Electro Medical Systems SA  
Ch. de la Vuarpillière 31  
CH-1260 Nyon - SWITZERLAND

Distributed by:  
Boston Scientific Corporation  
300 Boston Scientific Way  
Marlborough, MA 01752  
www.bostonscientific.com

**Ordering Information**  
**1.888.272.1001**

© 2018 Boston Scientific Corporation  
or its affiliates. All rights reserved.

URO-520103-AA APR 2018

For more information, visit  
[BostonScientific.com/Trilogy](https://www.bostonscientific.com/Trilogy)

\* Optional feature. Check for availability.

Caution: U.S. Federal law restricts this device to sale by or on the order of a physician.

All trademarks are the property of their respective owners.